

Laser-Induced Breakdown Spectroscopy for the Detection and Diagnosis of Bacterial Pathogens in Blood, Urine, and Cerebrospinal Fluid

Steven J. Rehse, Emma Blanchette, Emily Tracey, Caroline Alionte, Hadia Malik, August Baughan, Isabella Arthur, Jasmine Saad, Rachel Chevalier, Nicholas Bolton, Matteo Pontoni, Lauren Dmytrow, Abdullah Mustafa, and Mila Vasquez

Department of Physics, University of Windsor, Windsor, Ontario, Canada



Historical Notes

LIBS has been used to detect/identify bacteria since 2003.

Detection of bacteria by time-resolved laser-induced breakdown spectroscopy

Stéphane Morel, Nicolas Leone, Philippe Adam, and Jacques Amouroux

6184 APPLIED OPTICS / Vol. 42, No. 30 / 20 October 2003

Laser-Induced Breakdown Spectroscopy Detection and Classification of Biological Aerosols*

JOHN D. HYBL, GREGG A. LITHGOW, and STEVEN G. BUCKLEY†
 MIT Lincoln Laboratory, 244 Wood Street, Lexington, Massachusetts 02420-9108 (J.D.H.); and Department of Mechanical and Aerospace Engineering, University of California, San Diego, California 92093 (G.A.L., S.G.B.)

Volume 57, Number 10, 2003




APPLIED SPECTROSCOPY 1207

Laser-induced breakdown spectroscopy of bacterial spores, molds, pollens, and protein: initial studies of discrimination potential


Alan C. Samuels, Frank C. DeLucia, Jr., Kevin L. McNesby, and Andrzej W. Miziolek

20 October 2003 / Vol. 42, No. 30 / APPLIED OPTICS 6205



Bacillus globigii BG-1

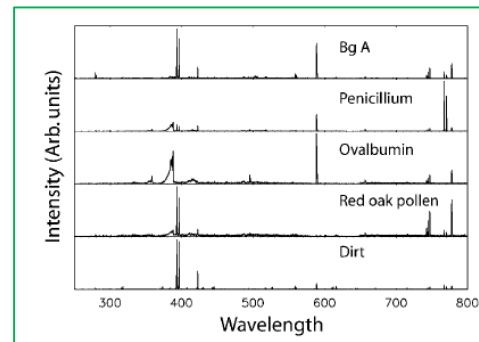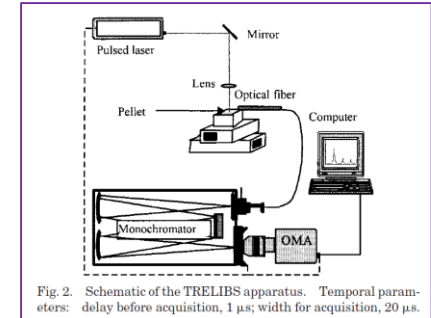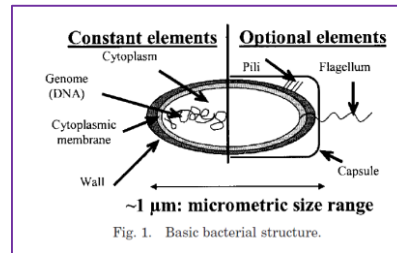

Bacillus subtilis var. niger

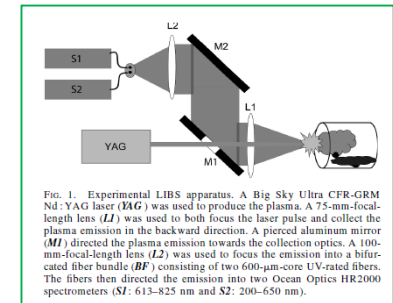

Bacillus thuringiensis subsp. *kurstaki*, *Bacillus subtilis* (also known as *Bacillus globigii*) subsp. *niger*, and *Bacillus cereus* 6E1

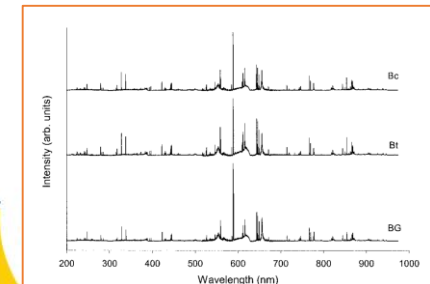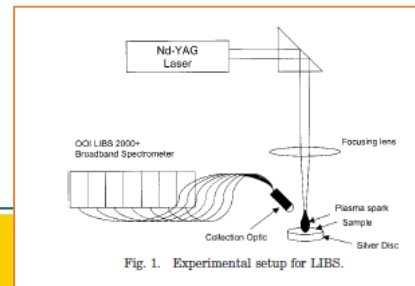

Historical Notes


Significant achievements followed quickly.

Anal. Chem. **2005**, *77*, 631–638

Feasibility of Detection and Identification of Individual Bioaerosols Using Laser-Induced Breakdown Spectroscopy

P. B. Dixon and D. W. Hahn*

Department of Mechanical and Aerospace Engineering, University of Florida, Gainesville, Florida 32611-6300

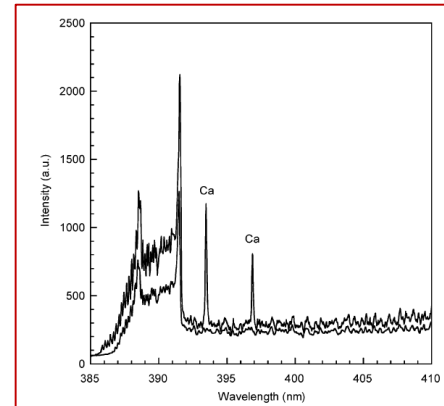

Detection of single *Bacillus* spores

Figure 3. Ensemble-averaged spectrum of 40 identified individual hits of *B. atrophaeus* spores, along with a representative ensemble-averaged spectrum corresponding to the absence of any spores. Spectra have the same intensity scale.

JOURNAL OF APPLIED PHYSICS **99**, 084701 (2006)


Femtosecond time-resolved laser-induced breakdown spectroscopy for detection and identification of bacteria: A comparison to the nanosecond regime

Matthieu Baudelet, Laurent Guyon, Jin Yu,^{a)} and Jean-Pierre Wolf
Laboratoire de Spectrométrie Ionique et Moléculaire, UMR CNRS 5579, Université Claude Bernard-Lyon 1, 43, Boulevard du 11 Novembre 1918, F-69622 Villeurbanne Cedex, France

Tanguy Amodeo and Emeric Fréjafon
Institut National de l'Environnement Industriel et des Risques (INERIS), Parc technologique ALATA, BP2, 60550 Verneuil-en-Halatte

Patrick Laloi
Laboratoire de Microbiologie et Génétique, UMR CNRS 5122, Université Claude Bernard-Lyon 1, 43, Boulevard du 11 Novembre 1918, F-69622 Villeurbanne Cedex, France

Escherichia coli and *Bacillus subtilis* with fs-LIBS

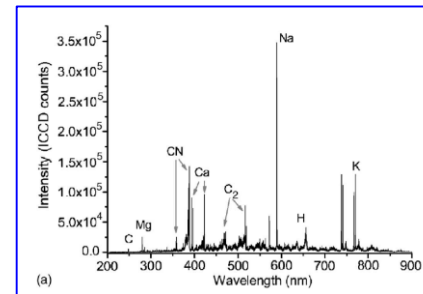

Discrimination of three strains of *Escherichia coli*, mold, and *Candida albicans* yeast


APPLIED PHYSICS LETTERS **90**, 163901 (2007)

Escherichia coli identification and strain discrimination using nanosecond laser-induced breakdown spectroscopy

Jonathan Diedrich and Steven J. Rehse^{a)}
Department of Physics and Astronomy, Wayne State University, Detroit, Michigan 48201

Sunil Palchaudhuri
Department of Immunology and Microbiology, Wayne State University, Detroit, Michigan 48201

(Received 26 January 2007; accepted 15 March 2007; published online 16 April 2007)

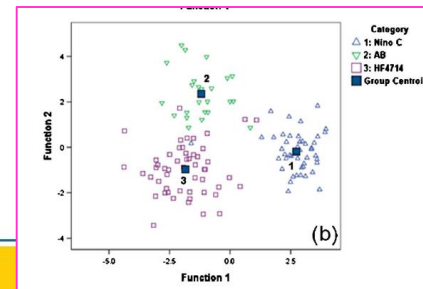

The Reasons?

There did not exist at the time, and there still is not now, a diagnostic test to detect and diagnose bacterial pathogens that can be:

- highly-flexible (multiple pathogens)
- very rapid (under 5 minutes)
- robust (field portable)
- relatively simple (minimal prep, performed by a non-expert)
- few to no consumables (biodegradation, cost)
- sensitive
- specific

Bioweapons / bioterrorism

Drinking water

Environmental surveillance

Food / beverage preparation

Medical diagnosis



To the present

Since then, progress has been extensive on:

- the number / variety of species
- sample preparation methods
- hardware setups
- chemometric algorithms
- number of bacteria detectable
- machine learning



To the present

Table 2

A summary of all the bacterial species/strains tested with LIBS.

Micro-organism	Reference	Form	Chemometric utilized	Laser wavelength ^a
<i>Acinetobacter baumannii</i> ATCC BAA-1789	[66]	Colony on blood agar	PCA/PLS1	1064
<i>Acinetobacter baylyi</i>	[48]	Pellet, freeze-dried powder	Hyperspace projection of trace elements	810 (fs)
<i>Acinetobacter calcoaceticus</i> [FJ816073] ^b	[71]	Colony on glass slide	PCA/PLS-RA	800 (fs)
<i>Arhodomonas</i> sp. [EU308280]	[71]	Colony on glass slide	PCA/PLS-RA	800 (fs)
<i>Bacillus anthracis</i> var. Sterne	[44]	Thin lawn on nylon filter	None	1064
<i>Bacillus anthracis</i> var. Sterne	[38]	Thin lawn on agar, glass slide	PCA/PLS1	1064
<i>Bacillus atrophaeus</i>	[41]	Spore, aerosol stream	None	1064
<i>Bacillus atrophaeus</i>	[52,37]	Dried film on Al disk, steel disk, polycarbonate disk	NN, MLSRA, PLS-DA	1064
<i>Bacillus atrophaeus</i>	[57]	Pellet, freeze-dried powder	SVM	1064
<i>Bacillus aureus</i>	[35]	Spore, EDB trap	None	355
<i>Bacillus cereus</i> 6E1	[29,32]	Thin lawn on silver membrane filter	PCA, linear correlation, and SIMCA	1064
<i>Bacillus cereus</i> ATCC 14603	[57]	Pellet, freeze-dried powder	SVM	1064
<i>Bacillus globigii</i> ^d BG-1	[39,40]	Pellet, freeze-dried powder	None	1064
<i>Bacillus globigii</i> BG-1	[40]	Spore, aerosol stream	None	1064
<i>Bacillus globigii</i> BG-2	[39,40]	Pellet, freeze-dried powder	None	1064
<i>Bacillus globigii</i> BG-2	[40]	Spore, aerosol stream	None	1064
<i>Bacillus globigii</i> var. niger	[29,42,32]	Thin lawn on silver membrane filter	PCA, linear correlation, and SIMCA	1064
<i>Bacillus globigii</i> var. niger	[30]	Continually refreshed dense aerosol cloud (from powder) and aerosol stream	PCA	1064
<i>Bacillus globigii</i> var. niger	[35,68]	Powder on double-sided sticky tape	No, linear correlation, PCA, PLS-DA	1064 × 2 (DP)
<i>Bacillus globigii</i> 168	[43]	Colony (wet) on LB medium	None	532
<i>Bacillus globigii</i>	[47]	Thin film lawn on cellulose nitrate membrane filter	None	810 (fs), 1064
<i>Bacillus globigii</i>	[48]	Pellet, freeze-dried powder	Hyperspace projection of trace elements	810 (fs)
<i>Bacillus globigii</i>	[31]	Dried powder on solid substrate	PCA, HCA, PCA + LDA	1064
<i>Bacillus globigii</i> ATCC 23857	[66]	Colony on blood agar	PCA/PLS1	1064
<i>Bacillus megaterium</i> QM B1551	[43]	Colony (wet) on LB medium	None	532
<i>Bacillus megaterium</i> PV361	[43]	Colony (wet) on LB medium	None	532
<i>Bacillus stearothersophilus</i> ATCC 12979	[57]	Pellet, freeze-dried powder	SVM	1064
<i>Bacillus thurengensis</i>	[39,40]	Pellet, freeze-dried powder	None	1064
<i>Bacillus thurengensis</i> var. kurstaki	[29,32]	Thin lawn on silver membrane filter	PCA, linear correlation, and SIMCA	1064
<i>Bacillus thurengensis</i> var. kurstaki	[44]	Thin lawn on nylon filter	None	1064
<i>Bacillus thurengensis</i> T34	[43]	Colony (wet) on LB medium	None	532
<i>Bacillus thuringiensis</i> ATCC 51912	[57]	Pellet, freeze-dried powder	SVM	1064
<i>Bacillus</i> sp. [GQ392044]	[71]	Colony on glass slide	PCA/PLS-RA	800 (fs)
<i>Bacillus</i> sp. [GQ226038]	[71]	Colony on glass slide	PCA/PLS-RA	800 (fs)
<i>Bacillus</i> sp. [HM026606]	[71]	Colony on glass slide	PCA/PLS-RA	800 (fs)
<i>Enterobacter cloacae</i> [FJ194527]	[71]	Colony on glass slide	PCA/PLS-RA	800 (fs)
<i>Enterobacter cloacae</i> ATCC 13047	[64,67]	Thin lawn on nutrient-free agar	DFA, PLS-DA	1064
<i>Enterobacter</i> sp. [CP000653]	[71]	Colony on glass slide	PCA/PLS-RA	800 (fs)
<i>Enterobacter</i> sp. [GU586319]	[71]	Colony on glass slide	PCA/PLS-RA	800 (fs)

Specimen presentation

In my research, I am almost completely focused on developing a realistic clinical testing methodology.

Current methods of LIBS bacterial analysis are completely inadequate for this purpose.

Most samples are obtained as fluid specimens, need to perform LIBS on that to be truly rapid.



Outline for rest of talk

1. Explain how we perform LIBS on bacteria in my lab in general

2. Three case studies

2a. Using the same sample preparation and data analysis


- **Blood**

- **Urine**

2b. Using a modified sample preparation and data analysis

- **Cerebrospinal fluid**



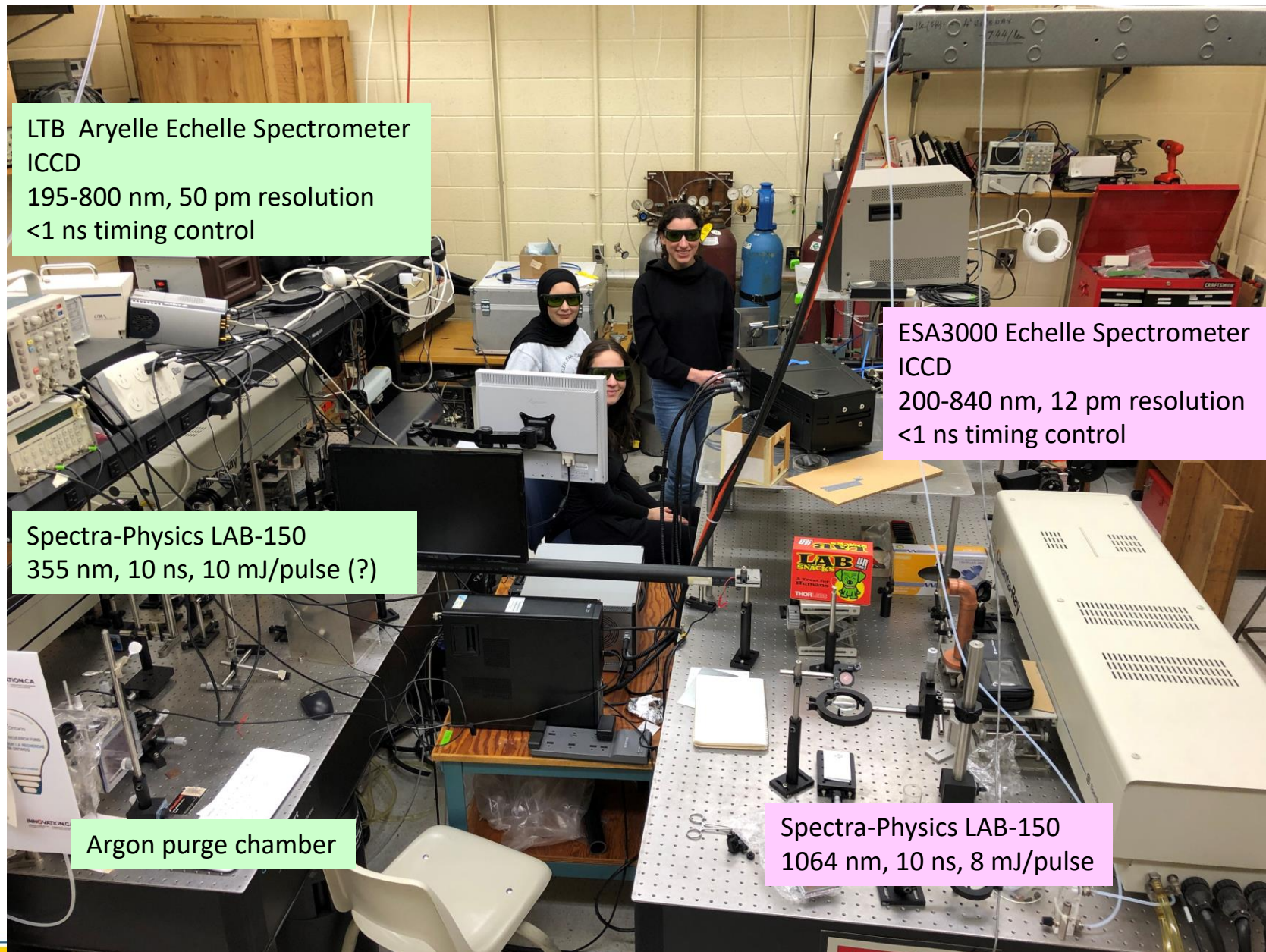

LTB Aryelle Echelle Spectrometer
ICCD
195-800 nm, 50 pm resolution
<1 ns timing control


ESA3000 Echelle Spectrometer
ICCD
200-840 nm, 12 pm resolution
<1 ns timing control

Spectra-Physics LAB-150
355 nm, 10 ns, 10 mJ/pulse (?)

Argon purge chamber

Spectra-Physics LAB-150
1064 nm, 10 ns, 8 mJ/pulse



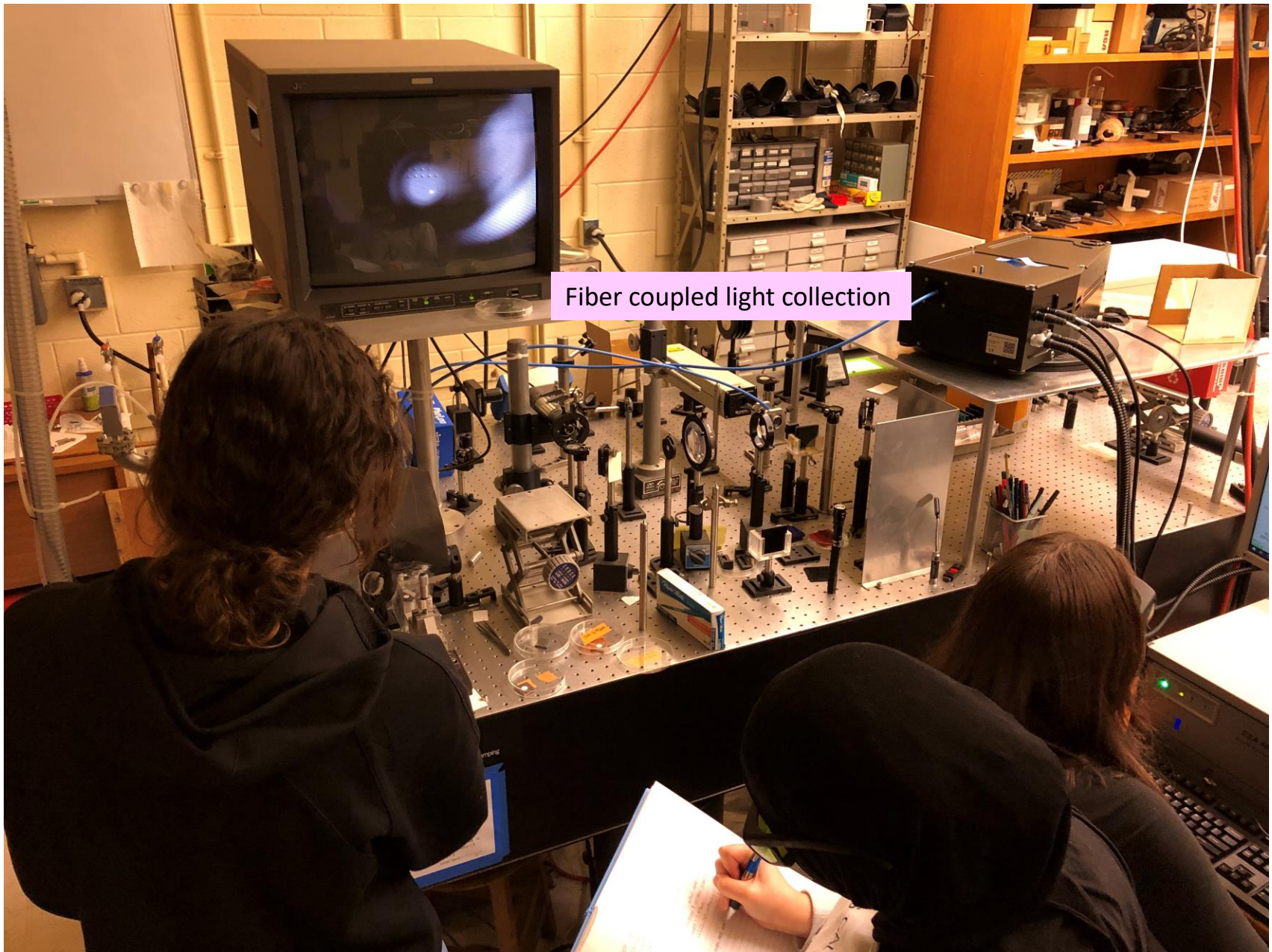


Fiber coupled light collection



Methodology – Bacterial Growth & Sample Prep




Bacteria cultured on agar plate

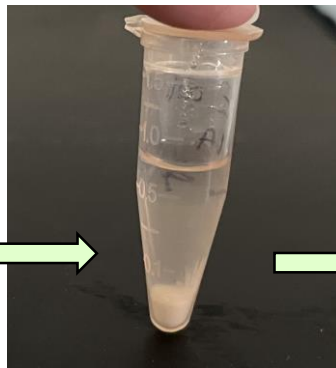


Harvested bacteria stored in MQH₂O

Bacteria-negative specimens collected from hospital

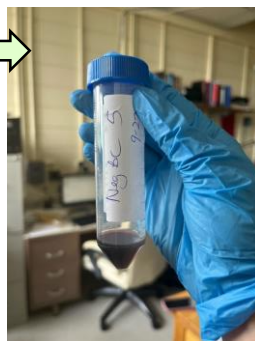


100 μ L fluid specimen + 100 μ L bacterial suspension pipetted into reassembled apparatus and sealed in centrifuge tube

9.5 mm diameter disposable filter

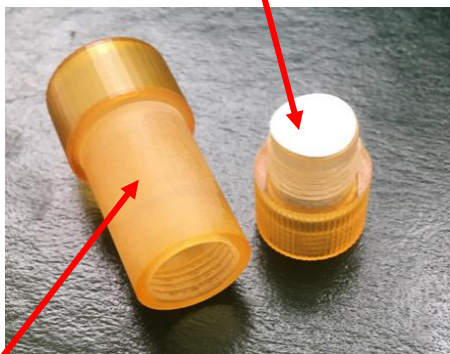

Disassembled 3D printed centrifuge insert

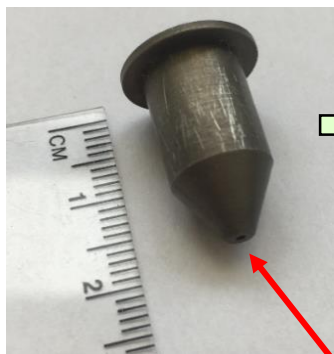

Aluminum cone with 1 mm diameter aperture at apex



Centrifuge insert assembled with filter & cone



Methodology – Bacterial Growth & Sample Prep

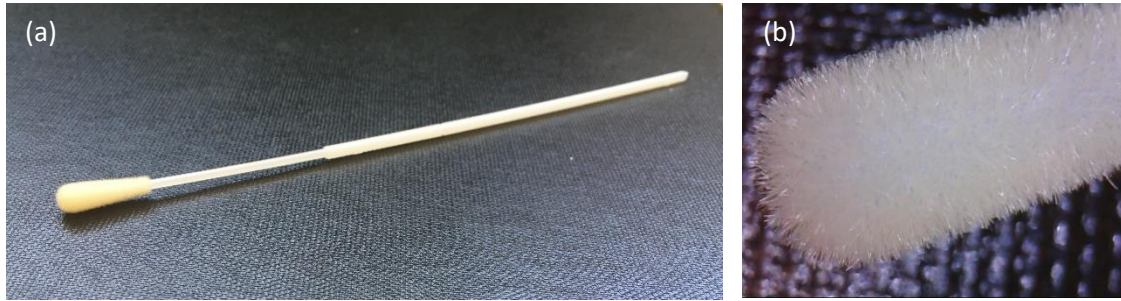


Figure 7.1: (a) Flocked swab used in this work. (b) Flocked swab zoomed-in on the tip.


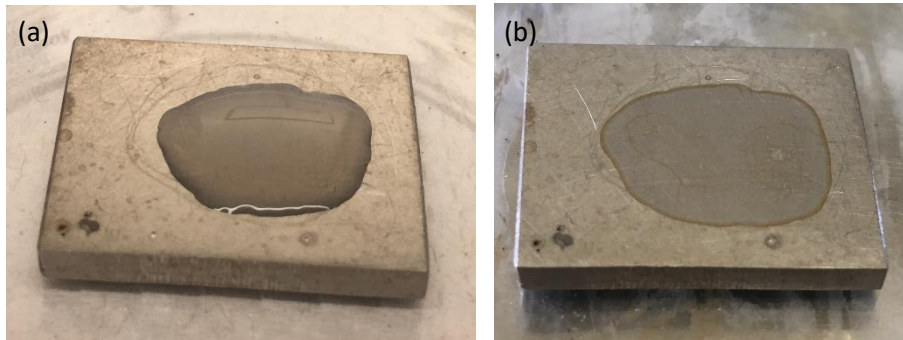


Figure 7.5: (a) 100 μL of *E. coli* pipetted onto surface of metal plate. (b) Metal plate after heated on hot-plate for 2 minutes 20 seconds at 200 $^{\circ}\text{C}$. Water has evaporated and film of bacteria is observed.




Flocked swab was pre-wet with 10 μL of deionized and swabbed over surface

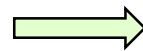

The swab was placed in 1 mL of deionized water and vortexed for 15 seconds. The entire 1 mL was deposited on a filter paper using the insert and metal cone.




Methodology – Bacterial Growth & Sample Prep

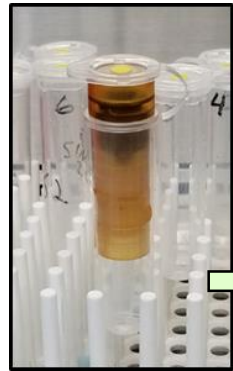


Assembled centrifuge tube insert


Samples centrifuged at 5000 RPM for five minutes

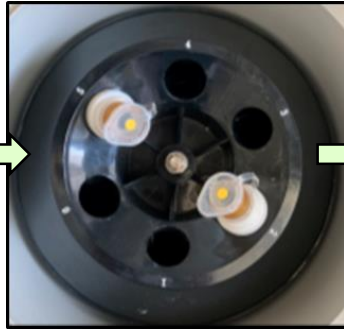



Insert disassembled, filter mounted on steel with double-sided sticky tape

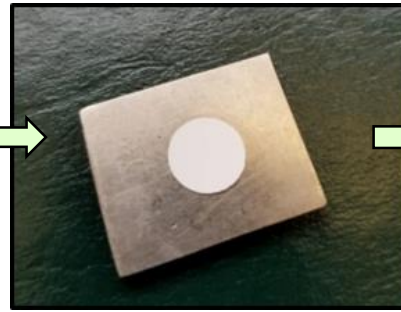

Sample ablated in Ar environment

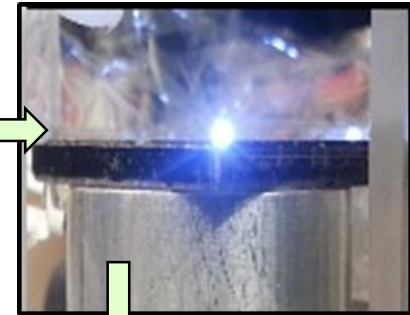

Filter after centrifugation; bacteria concentrated at the center. We obtain 25-30 LIBS spectra per deposition (could be more)

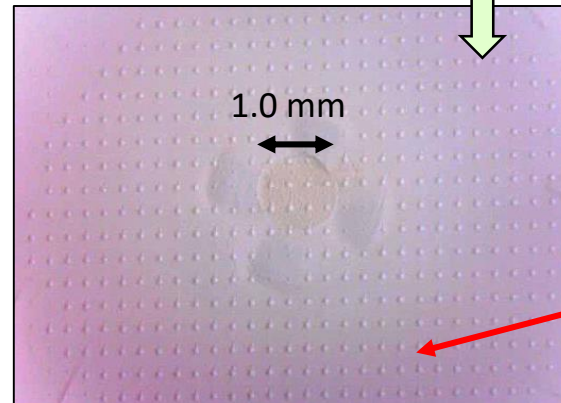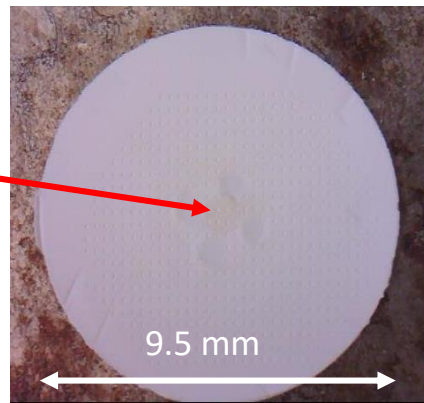


Magnification of the filter with a rastered array of laser ablation sites

Images of filter after data acquisition




Methodology – Bacterial Growth & Sample Prep

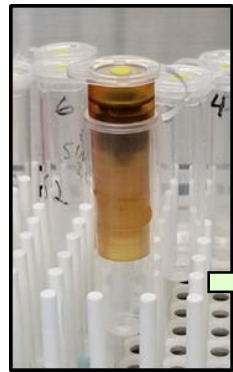


Assembled centrifuge tube insert


Samples centrifuged at 5000 RPM for five minutes

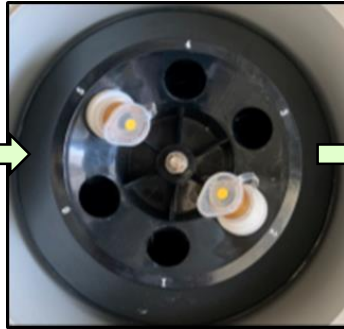



Insert disassembled, filter mounted on steel with double-sided sticky tape

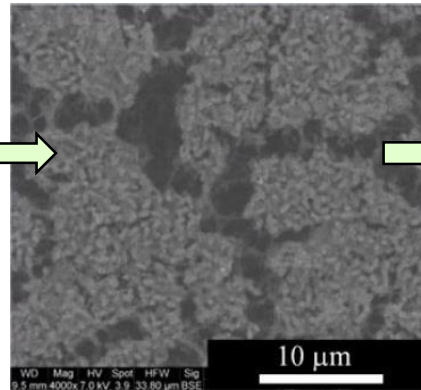

Sample ablated in Ar environment

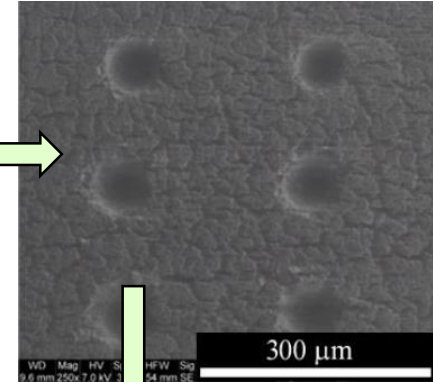

Filter after centrifugation; bacteria concentrated at the center. We obtain 25-30 LIBS spectra per deposition (could be more)

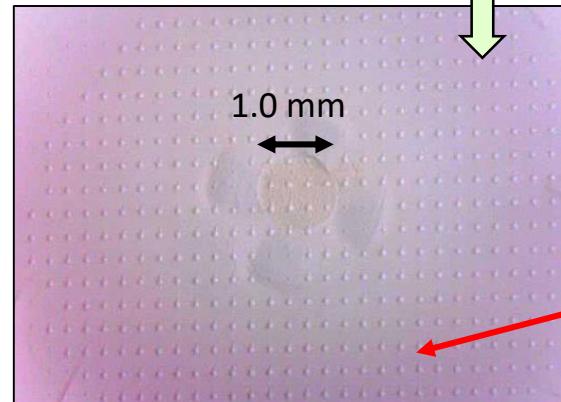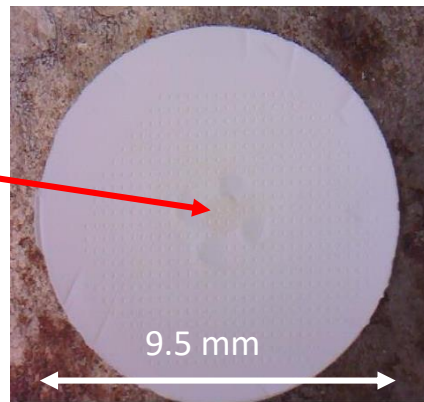

Magnification of the filter with a rastered array of laser ablation sites

Images of filter after data acquisition




Methodology – Bacterial Growth & Sample Prep




Samples centrifuged at

Insert disassembled, filter mounted on steel with

Sample ablated in

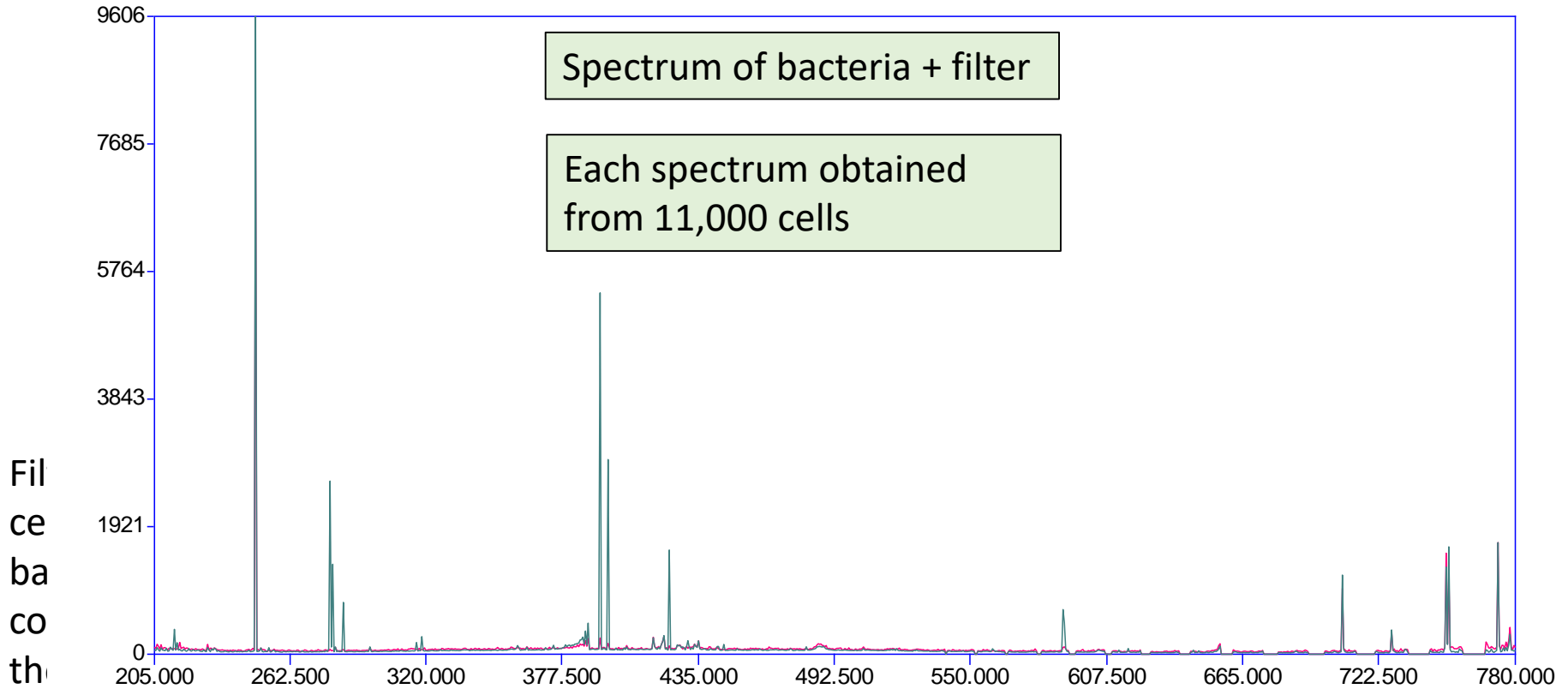

obtain 25-30 LIBS spectra per deposition (could be more)

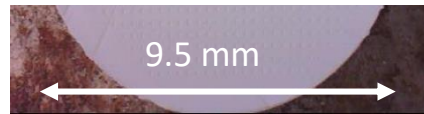


Images of filter after data acquisition



Spectrum From Bacteria

An elemental assay of the bacterial cell composition!

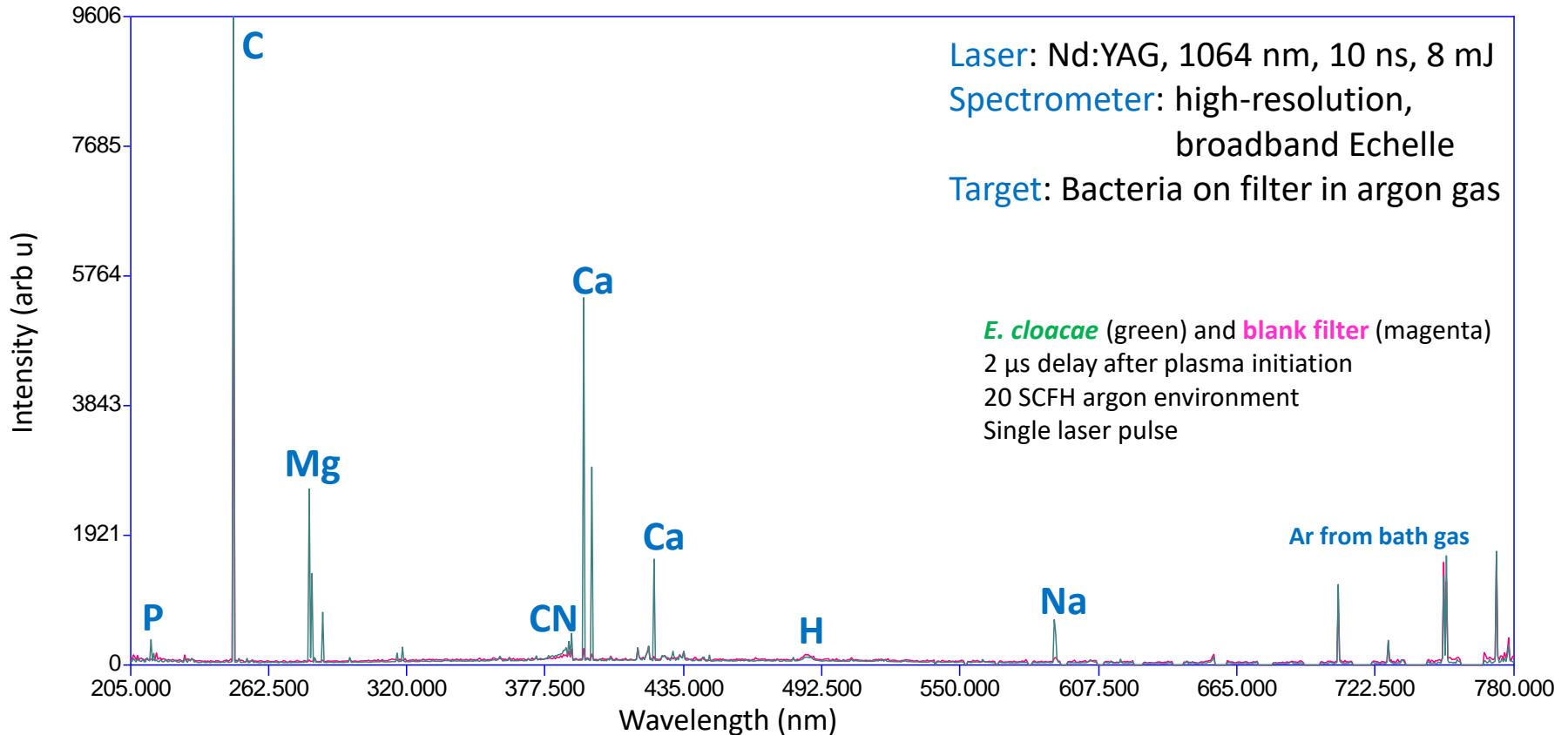

Suggests a real-time method for pathogenic bacterial diagnosis.



Detecting Species in Fluids with Chemometrics


“Ratio Model” consists of:

- 15 emission line intensities from P, C, Mg, Ca, and Na and 92 simple ratios
- 107 total independent variables

PLS-DA on
Ratio Model

PLS-DA implemented with PLS_Toolbox (Eigenvector, Inc)

- Run on desktop PC, takes seconds
- Always perform a two-class test and remove entire filters at a time for testing

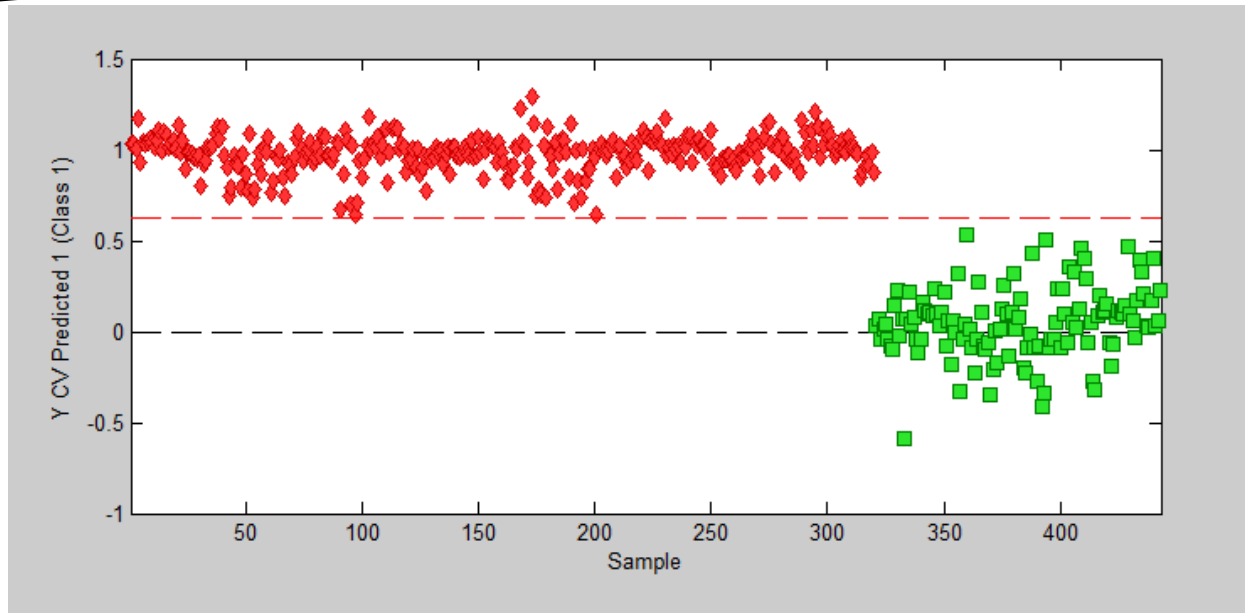

Diagnosing Species in Fluids with Machine Learning

“Ratio Model” consists of:

- 15 emission line intensities and 92 simple ratios
- 107 independent variables

DFA on
Ratio Model

ANN on Ratio
Model

PCA-ANN on Ratio
Model

PCA-ANN on Full
Spectrum Data*

- Started using the “whole spectrum” from 200 nm - 590 nm. 42,000 variables.
- Perform unsupervised PCA first (implemented in Python), reduce to 10 PC's.
- Models are trained on 80% of single shot data, 20% reserved for testing. (~15 seconds).


PCA-ANN With Full Spectrum Data			
	E. coli	S. aureus	E. cloacae
Sensitivity	98.04 %	93.27 %	91.23 %
Specificity	97.71 %	97.22 %	96.12 %
Classification Error	2.13 %	4.28 %	6.33 %

*The full spectrum spans 200 nm – 760 nm, but no lines of interest > 590 nm



Diagnosing Species in Water with Machine Learning

- PCA-ANN custom written in Python
- 42,000 independent variables \Rightarrow 10 principal components (new independent variables)
- Currently only one hidden layer
- ANN algorithm implemented to optimize the patience and # of hidden nodes
- Models are trained on 80% of single shot data, 20% reserved for testing.
- After PCA, spectra were randomized to demonstrate the ANN was not fitting noise (random classification of cases)

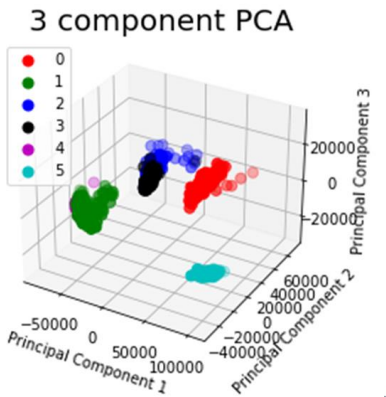


Python (Python 3.8)

PCA-ANN With Full Spectrum Data

	E. coli	S. aureus	E. cloacae
Sensitivity	98.04 %	93.27 %	91.23 %
Specificity	97.71 %	97.22 %	96.12 %
Classification Error	2.13 %	4.28 %	6.33 %

```

100: # PCA-ANN With Full Spectrum Data
101: # PCA-ANN With Full Spectrum Data
102: # PCA-ANN With Full Spectrum Data
103: # PCA-ANN With Full Spectrum Data
104: # PCA-ANN With Full Spectrum Data
105: # PCA-ANN With Full Spectrum Data
106: # PCA-ANN With Full Spectrum Data
107: # PCA-ANN With Full Spectrum Data
108: # PCA-ANN With Full Spectrum Data
109: # PCA-ANN With Full Spectrum Data
110: # PCA-ANN With Full Spectrum Data
111: # PCA-ANN With Full Spectrum Data
112: # PCA-ANN With Full Spectrum Data
113: # PCA-ANN With Full Spectrum Data
114: # PCA-ANN With Full Spectrum Data
115: # PCA-ANN With Full Spectrum Data
116: # PCA-ANN With Full Spectrum Data
117: # PCA-ANN With Full Spectrum Data
118: # PCA-ANN With Full Spectrum Data
119: # PCA-ANN With Full Spectrum Data
120: # PCA-ANN With Full Spectrum Data
121: # PCA-ANN With Full Spectrum Data
122: # PCA-ANN With Full Spectrum Data
123: # PCA-ANN With Full Spectrum Data
124: # PCA-ANN With Full Spectrum Data
125: # PCA-ANN With Full Spectrum Data
126: # PCA-ANN With Full Spectrum Data
127: # PCA-ANN With Full Spectrum Data
128: # PCA-ANN With Full Spectrum Data
129: # PCA-ANN With Full Spectrum Data
130: # PCA-ANN With Full Spectrum Data
131: # PCA-ANN With Full Spectrum Data
132: # PCA-ANN With Full Spectrum Data
133: # PCA-ANN With Full Spectrum Data
134: # PCA-ANN With Full Spectrum Data
135: # PCA-ANN With Full Spectrum Data
136: # PCA-ANN With Full Spectrum Data
137: # PCA-ANN With Full Spectrum Data
138: # PCA-ANN With Full Spectrum Data
139: # PCA-ANN With Full Spectrum Data
140: # PCA-ANN With Full Spectrum Data
141: # PCA-ANN With Full Spectrum Data
142: # PCA-ANN With Full Spectrum Data
143: # PCA-ANN With Full Spectrum Data
144: # PCA-ANN With Full Spectrum Data
145: # PCA-ANN With Full Spectrum Data
146: # PCA-ANN With Full Spectrum Data
147: # PCA-ANN With Full Spectrum Data
148: # PCA-ANN With Full Spectrum Data
149: # PCA-ANN With Full Spectrum Data
150: # PCA-ANN With Full Spectrum Data
151: # PCA-ANN With Full Spectrum Data
152: # PCA-ANN With Full Spectrum Data
153: # PCA-ANN With Full Spectrum Data
154: # PCA-ANN With Full Spectrum Data
155: # PCA-ANN With Full Spectrum Data
156: # PCA-ANN With Full Spectrum Data
157: # PCA-ANN With Full Spectrum Data
158: # PCA-ANN With Full Spectrum Data
159: # PCA-ANN With Full Spectrum Data
160: # PCA-ANN With Full Spectrum Data
161: # PCA-ANN With Full Spectrum Data
162: # PCA-ANN With Full Spectrum Data
163: # PCA-ANN With Full Spectrum Data
164: # PCA-ANN With Full Spectrum Data
165: # PCA-ANN With Full Spectrum Data
166: # PCA-ANN With Full Spectrum Data
167: # PCA-ANN With Full Spectrum Data
168: # PCA-ANN With Full Spectrum Data
169: # PCA-ANN With Full Spectrum Data
170: # PCA-ANN With Full Spectrum Data
171: # PCA-ANN With Full Spectrum Data
172: # PCA-ANN With Full Spectrum Data
173: # PCA-ANN With Full Spectrum Data
174: # PCA-ANN With Full Spectrum Data
175: # PCA-ANN With Full Spectrum Data
176: # PCA-ANN With Full Spectrum Data
177: # PCA-ANN With Full Spectrum Data
178: # PCA-ANN With Full Spectrum Data
179: # PCA-ANN With Full Spectrum Data
180: # PCA-ANN With Full Spectrum Data
181: # PCA-ANN With Full Spectrum Data
182: # PCA-ANN With Full Spectrum Data
183: # PCA-ANN With Full Spectrum Data
184: # PCA-ANN With Full Spectrum Data
185: # PCA-ANN With Full Spectrum Data
186: # PCA-ANN With Full Spectrum Data
187: # PCA-ANN With Full Spectrum Data
188: # PCA-ANN With Full Spectrum Data
189: # PCA-ANN With Full Spectrum Data
190: # PCA-ANN With Full Spectrum Data
191: # PCA-ANN With Full Spectrum Data
192: # PCA-ANN With Full Spectrum Data
193: # PCA-ANN With Full Spectrum Data
194: # PCA-ANN With Full Spectrum Data
195: # PCA-ANN With Full Spectrum Data
196: # PCA-ANN With Full Spectrum Data
197: # PCA-ANN With Full Spectrum Data
198: # PCA-ANN With Full Spectrum Data
199: # PCA-ANN With Full Spectrum Data
200: # PCA-ANN With Full Spectrum Data
    
```

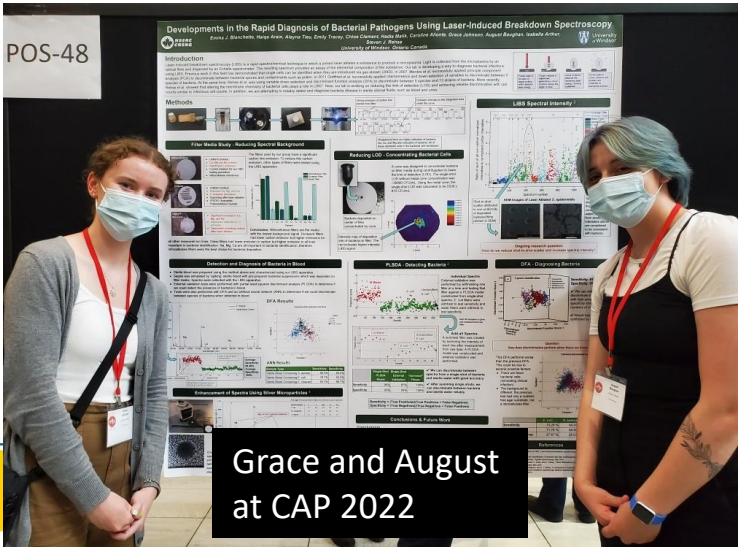

Blood (sepsis)

Sepsis (bacterial blood infection) is a highly dangerous infection and is a common nosocomial (hospital acquired) infection.

Sepsis is not always caused by bacterial infection, but many cases are bacterial in origin.

A blood culture is required for the diagnosis of **sepsis**, but not all blood cultures will produce a positive test result.

Antibiotics are the preferred treatment for **sepsis**, but even with the best treatment the mortality for patients who have reached septic shock is no better than 50%.

The time it takes to initiate effective antimicrobial therapy is the single strongest predictor of patient outcome, with every hour of delay increasing patient mortality.

E. Whitnack, Sepsis, in: N.C. Engleberg, V.J. DiRita, T. Dermody, M. Schaechter (Eds.), Mechanisms of Microbial Disease, 3rd ed, Lippincott Williams & Wilkins, 2007, pp. 564–572.

A. Kumar, D. Roberts, K.E. Wood, B. Light, J.E. Parrillo, et al., Duration of hypotension before initiation of effective antimicrobial therapy is the critical determinant of survival in human septic shock, Crit. Care Med. 34 (2006) 1589–1596.



Blood (sepsis)

The standard diagnostic blood cultures are slow and labor intensive, usually requiring a two-step procedure:

- 1) an initial culture taking up to 120 h
- 2) a 24-48 h subculture to identify the pathogens detected in a positive initial culture.



Automated blood culture system



Culturing and diagnosis done manually

J.C. Lagier, S. Edouard, I. Pagnier, O. Mediannikov, M. Drancourt, D. Raoult, Current and past strategies for bacterial culture in clinical microbiology, *Clin. Microbiol. Rev.* 28 (2015) 208–236, <https://doi.org/10.1128/cmr.00110-14>

D.J. Shin, N. Andini, K. Hsieh, S. Yang, T.-H. Wang, Emerging analytical techniques for rapid pathogen identification and susceptibility testing, *Ann. Rev. Anal. Chem.* 12 (2019) 41-67, <https://doi.org/10.1146/annurev-anchem-061318-115529>



Detecting Infection in Blood with a PLS-DA



Bacterial-negative specimens obtained from Windsor Regional Hospital with the anticoagulant sodium polyanetholesulfonate (SPS) present.

Five different blood samples from different patients were provided to account for the difference between patient physiology.


Specimens were vortex-mixed prior to handling to insure the redistribution of any suspended particles.

DETECTION IN BLOOD (individual spectra)

96.3% sensitivity – 600 bacterial spectra
 98.6% specificity – 206 blood spectra
 11,000 cells / spectrum

DETECTION IN BLOOD (spectra added)


19 of 19 filters positive, 100% sensitivity
 7 of 7 filters negative, 100% specificity

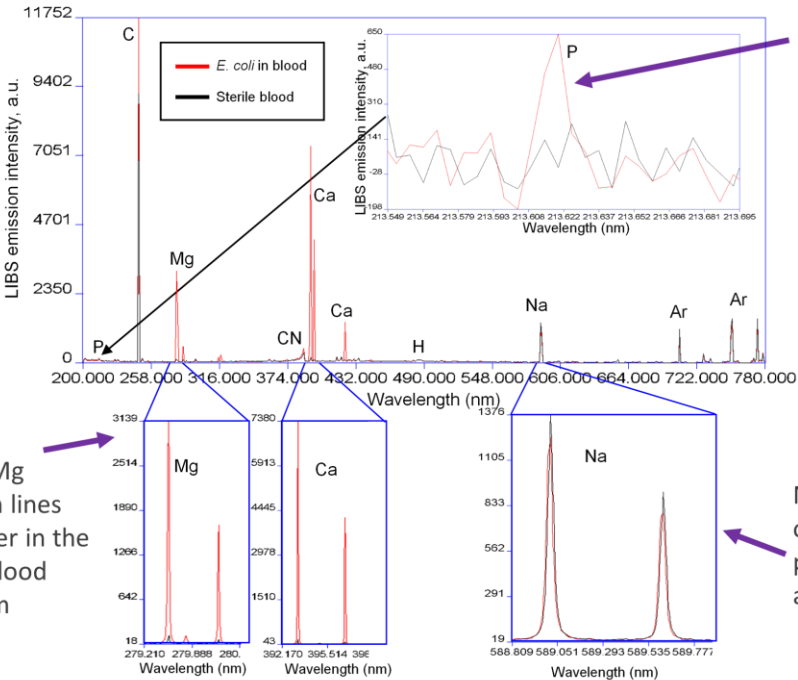

P line is small in the spiked blood, but not present in the sterile blood

Ca and Mg emission lines are higher in the spiked blood spectrum

Na emission observed is partially due to anticoagulant

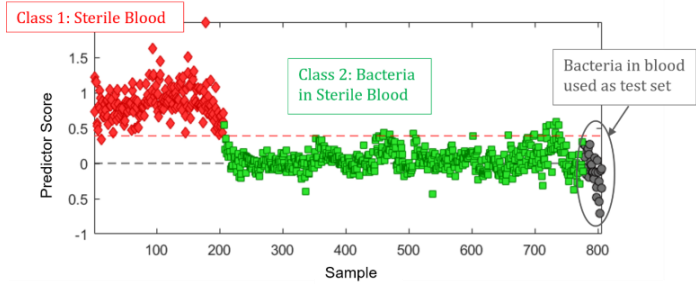

Prepared Samples:

- 5 filters each of *S. aureus*, *E. coli*, *E. cloacae*
- 4 filters of *P. aeruginosa*
- 7 filters of sterile blood
- 30 LIBS laser shots from each deposition



Diagnosing Infection in Blood with PCA-ANN




Bacterial-negative specimens obtained from Windsor Regional Hospital with the anticoagulant sodium polyanetholesulfonate (SPS) present.

Five different blood samples from different patients were provided to account for the difference between patient physiology.

Specimens were vortex-mixed prior to handling to insure the redistribution of any suspended particles.

PCA-ANN With Full Spectrum Data*				
	<i>S. aureus</i>	<i>E. coli</i>	<i>E. cloacae</i>	<i>P. aeruginosa</i>
Sensitivity	100 %	100 %	100 %	100 %
Specificity	100 %	100 %	100 %	100 %
Classification Error	0.00 %	0.00 %	0.00 %	0.00 %

*classification using 80:20 split

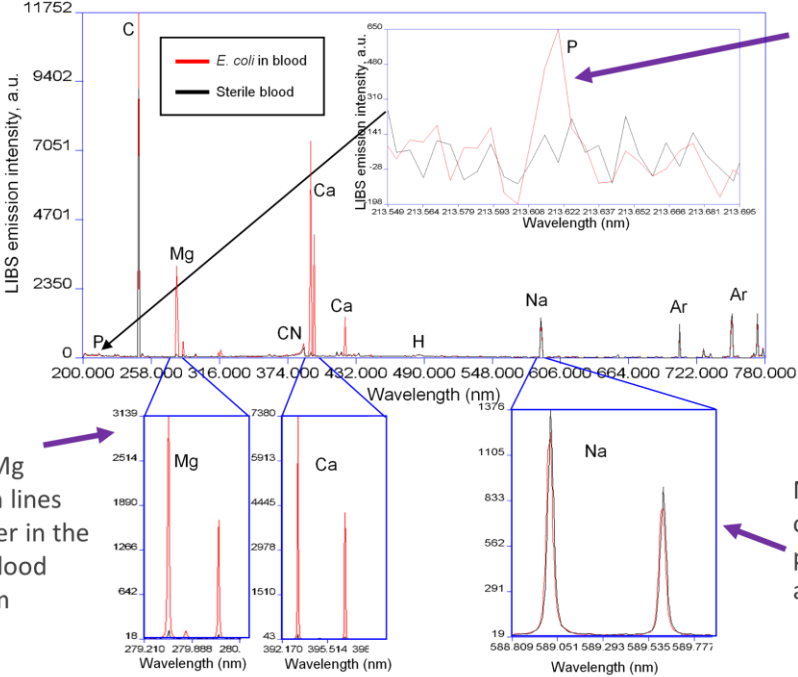

P line is small in the spiked blood, but not present in the sterile blood

Ca and Mg emission lines are higher in the spiked blood spectrum

Na emission observed is partially due to anticoagulant

Prepared Samples:

- 5 filters each of *S. aureus*, *E. coli*, *E. cloacae*
- 4 filters of *P. aeruginosa*
- 7 filters of sterile blood
- 30 LIBS laser shots from each deposition



Urine (UTI)

UTI (urinary tract infections) are one of the most common infections in adult women, with up to 50% of women having experienced at least one **UTI** in their lifetime and upwards of 10% of women experiencing at least one **UTI** annually.

In 2007, **UTI** resulted in approximately 8.6 million health care visits in the United States with an estimated cost of 1.6 billion US dollars.

In addition to the prevalence in the population, **UTI** are also the most common nosocomial infection, with up to 80% of hospital-acquired **UTI** associated with the use of a bladder catheter.



Even a single catheterization can lead to **UTI** due to contamination of the catheter tip.

T.K. Price, E.E. Hilt, T.J. Dune, E.R. Mueller, A.J. Wolfe, L. Brubaker, Urine trouble: should we think differently about UTI?, Int. Urogynecol. J. 29 (2018) 205-210,
S.M. Schappert, E.A. Rechtsteiner, Ambulatory medical care utilization estimates for 2007, Vital Health Stat. 13 169 (2011) 1-38, PMID: 21614897.

J.M.T. Barfor, K. Anson, Y.Hu, A.R.M. Coates, A model of catheter-associated urinary tract infection initiated by bacterial contamination of the catheter tip, BJU Int. 102 (2008) 67-74.



Urine (UTI)

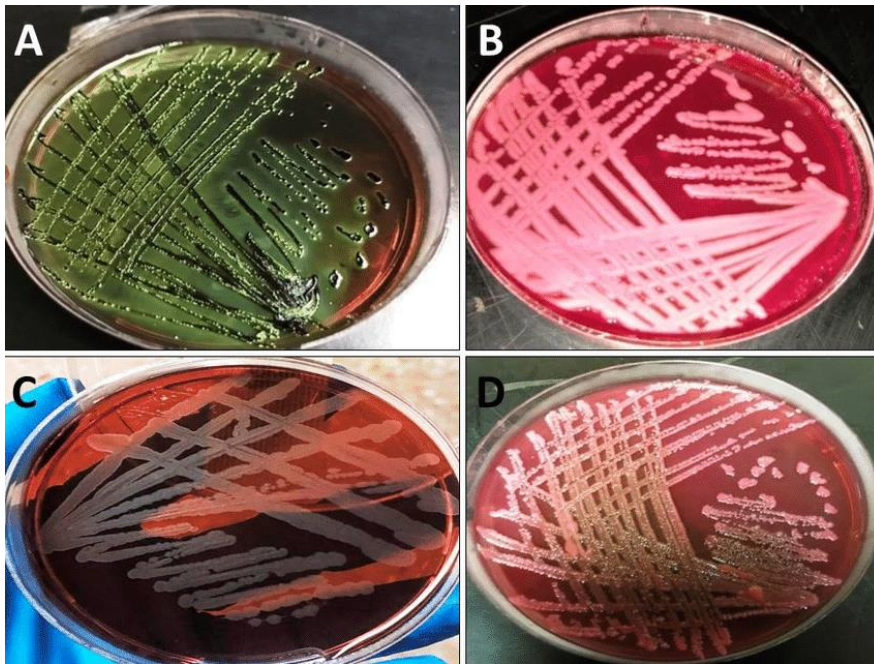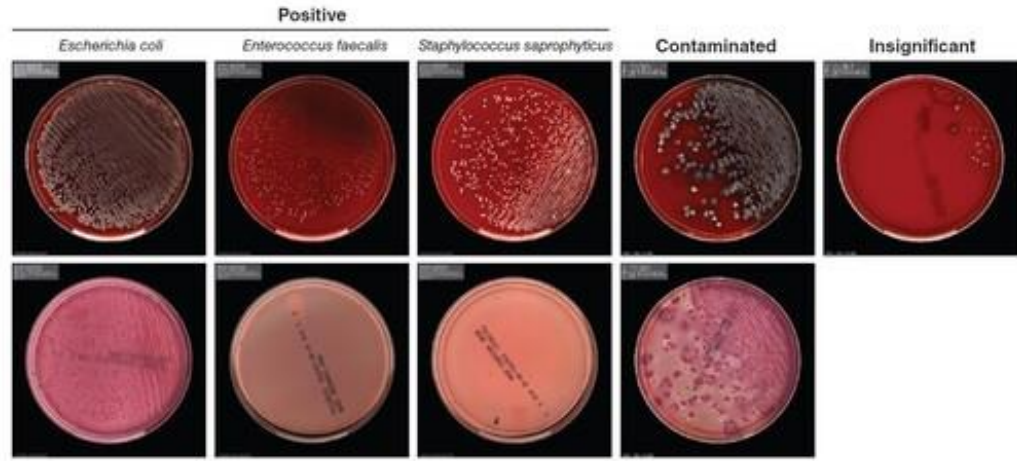

Detecting Infection in Urine with a PLS-DA



Bacterial-negative urine specimens obtained from Windsor Regional Hospital.

Urine specimens from four different patients were characterized to account for differences in patient physiology.


Specimens were vortex-mixed prior to handling to insure the redistribution of any suspended particles.

DETECTION IN URINE (individual spectra)

- 98.9% sensitivity – 360 bacterial spectra
- 100% specificity – 240 urine spectra
- 11,000 cells per spectrum

DETECTION IN URINE (spectra added)

- 12 of 12 filters positive, 100% sensitivity
- 8 of 8 filters negative, 100% specificity

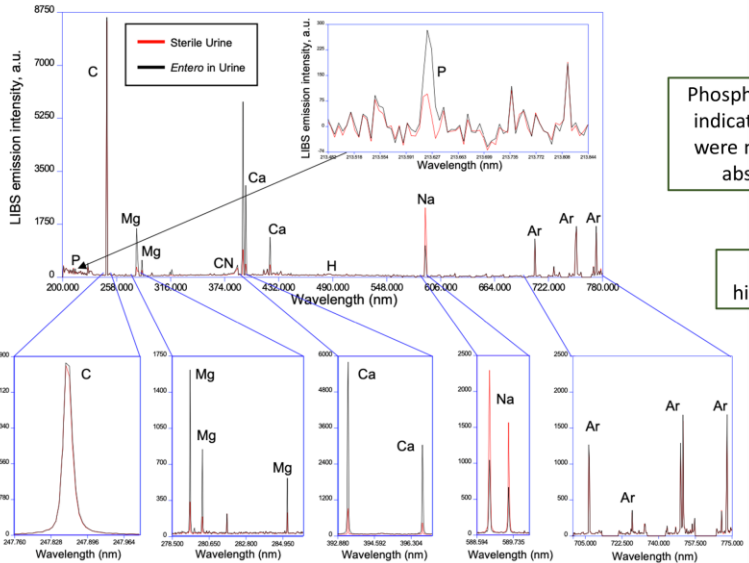

Strong carbon emission is predominantly from the nitrocellulose filter


Magnesium lines are greater in spiked urine, an indication of bacteria

Calcium lines higher in spiked urine samples, likely due to bacteria

Phosphorus lines are highly indicative of bacteria, and were never present in the absence of bacteria

Sodium emission higher in sterile urine

Argon lines are from the chamber environment where the samples are shot in

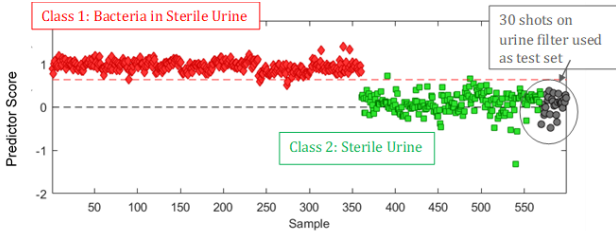

- Prepared Samples:**
- 4 filters each of *S. aureus*, *E. coli*, *E. cloacae*
 - 8 filters of sterile urine
 - 30 LIBS laser shots from each deposition



Diagnosing Infection in Urine with PCA-ANN




Bacterial-negative urine specimens obtained from Windsor Regional Hospital.

Urine specimens from four different patients were characterized to account for differences in patient physiology.

Specimens were vortex-mixed prior to handling to insure the redistribution of any suspended particles.

PCA-ANN With Full Spectrum Data*			
	<i>S. aureus</i>	<i>E. coli</i>	<i>E. cloacae</i>
Sensitivity	100 %	100 %	91.67 %
Specificity	100 %	95.83 %	100 %
Classification Error	0.00 %	2.09 %	4.17 %

*classification using 80:20 split

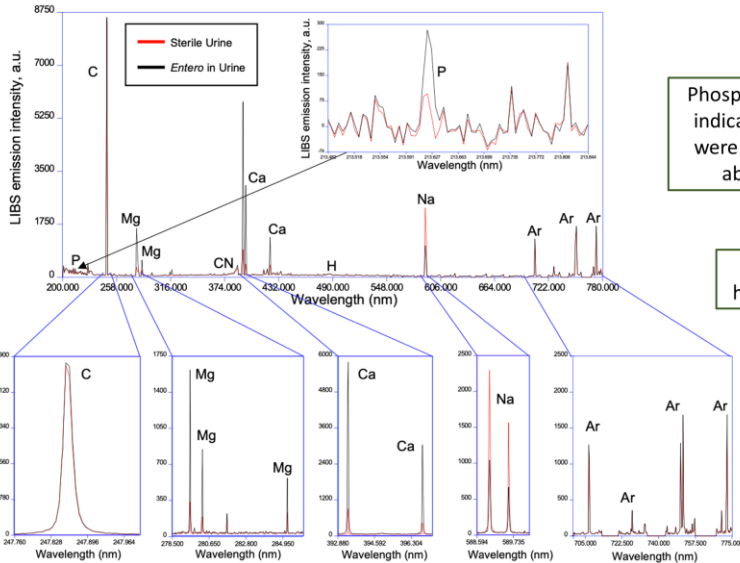

Strong carbon emission is predominantly from the nitrocellulose filter

Magnesium lines are greater in spiked urine, an indication of bacteria

Calcium lines higher in spiked urine samples, likely due to bacteria

Phosphorus lines are highly indicative of bacteria, and were never present in the absence of bacteria

Sodium emission higher in sterile urine

Argon lines are from the chamber environment where the samples are shot in

- Prepared Samples:**
- 4 filters each of *S. aureus*, *E. coli*, *E. cloacae*
 - 8 filters of sterile urine
 - 30 LIBS laser shots from each deposition

External Validation of PCA-ANN in Blood and Urine

- External validation done in both urine and blood (whole filters withheld)
- ANN model run 10 times per excluded filter
- Improvements need to be made on quality of spectra, not ML

average sensitivity of 74.6%

Urine	<i>E. coli</i>	<i>S. aureus</i>	<i>E. cloacae</i>	
Average Sensitivity	63.58 %	86.92 %	62.33 %	
Blood	<i>E. coli</i>	<i>S. aureus</i>	<i>E. cloacae</i>	<i>P. aeruginosa</i>
Average Sensitivity	93.73 %	61.73 %	91.60 %	95.0 %

average sensitivity of 85.5%

Urine

Sensitivity of each E. coli filter removed from the model to be externally validated. ANN run 10 times per filter.

<i>E. coli</i>	Predicted			Sensitivity
	<i>S. aureus</i>	<i>E. coli</i>	<i>E. cloacae</i>	
Sample #				
1	30	250	20	0.833
2	1	138	161	0.460
3	11	153	136	0.510
4	72	222	6	0.74
Sum	114	763	323	0.636

Blood

Sensitivity results for S. aureus filters removed from the model to be externally validated. ANN run 10 times per filter.

<i>S. aureus</i>	Predicted				Sensitivity
	<i>S. aureus</i>	<i>E. coli</i>	<i>E. cloacae</i>	<i>P. aeruginosa</i>	
Sample #					
1	50	54	86	110	0.167
2	300	0	0	0	1.000
3	1	27	266	6	0.033
4	290	10	0	0	0.967
5	285	13	0	2	0.950
Sum	926	104	352	118	0.617



Cerebrospinal Fluid (meningitis)

Cerebrospinal fluid (CSF) is the fluid the brain, cranium, and spinal cord are bathed in.

Meningitis is an infection of the brain and spinal cord's protective layers (cells in CSF).

Meningitis may be viral or bacterial in nature.

Delays and inappropriate antibiotic use can lead to severe complications, including irreversible brain damage or death.

The gold standard for diagnosis is a spinal tap (lumbar puncture) to remove **cerebrospinal fluid (CSF)** followed by culture, serology, and analysis.



Slane, Valori H. and Chandrashekar G. Unakal. "Tuberculous Meningitis." StatPearls, StatPearls Publishing, 18 November 2022.

Telano, Lauren N. and Stephen Baker. "Physiology, Cerebral Spinal Fluid." StatPearls, StatPearls Publishing, 4 July 2023.



Cerebrospinal Fluid (meningitis)

Spinal taps are rarely performed in Windsor.

Artificial cerebrospinal fluid was used in this study (aCSF).

Artificial cerebrospinal spinal fluid (aCSF) closely matches the electrolyte concentrations of actual CSF and is prepared from analytical grade reagents and high purity water.

Ion	Concentrations (mM)
Na ⁺	150
K ⁺	3.0
Ca ²⁺	1.4
Mg ²⁺	0.8
p ³⁻	1.0
Cl ⁻	155




New Apparatus Design

Two components, with a built-in concentration cone.


- ✓ Concentrates bacteria in a 1 mm diameter circle
- ✓ The ridge on the aperture of β piece prevents bacterial escape, allowing only water and ions to diffuse
- ✓ Indent created on filter is more easily identifiable by the LIBS apparatus and seals bacteria

Stronger bacterial signal and more reproducible data

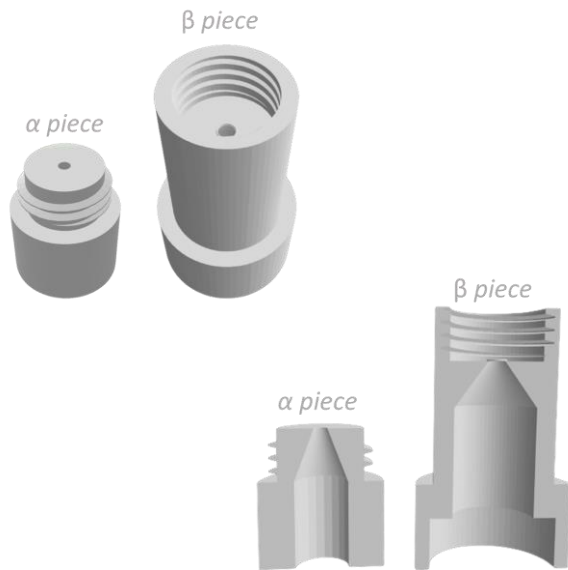

New Apparatus Data

*nitrocellulose filter
(0.45 μm pore size)*

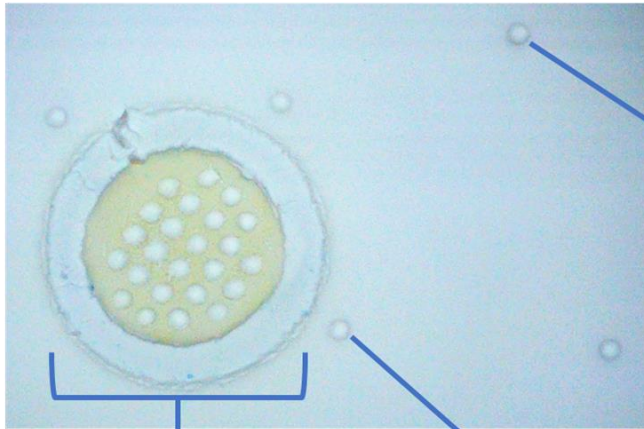

Laser shot
obtained far
from deposition
region and near
edge of filter or
“Far-Out”



Bacterial
deposition region
bounded by
custom cone
indent

Laser shot obtained
near to deposition
region or “Just-Out”




New Apparatus Data






Spectrum Obtained Inside Bacteria Deposition Region

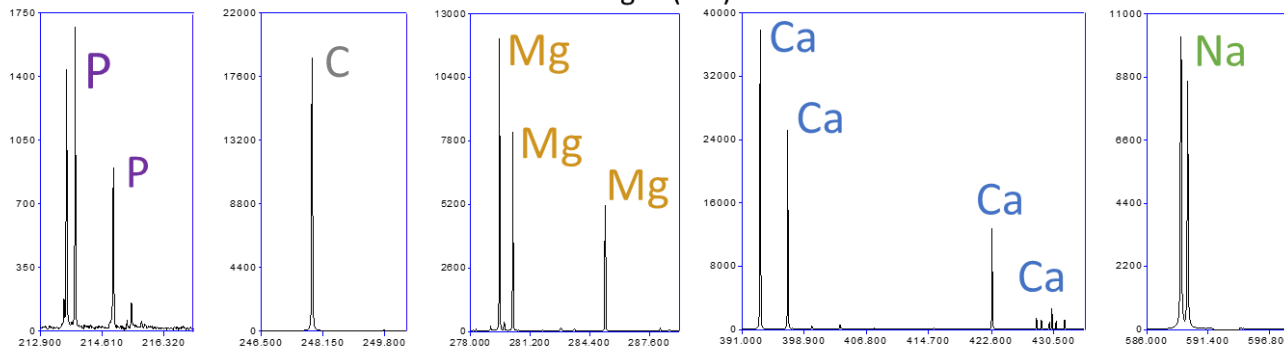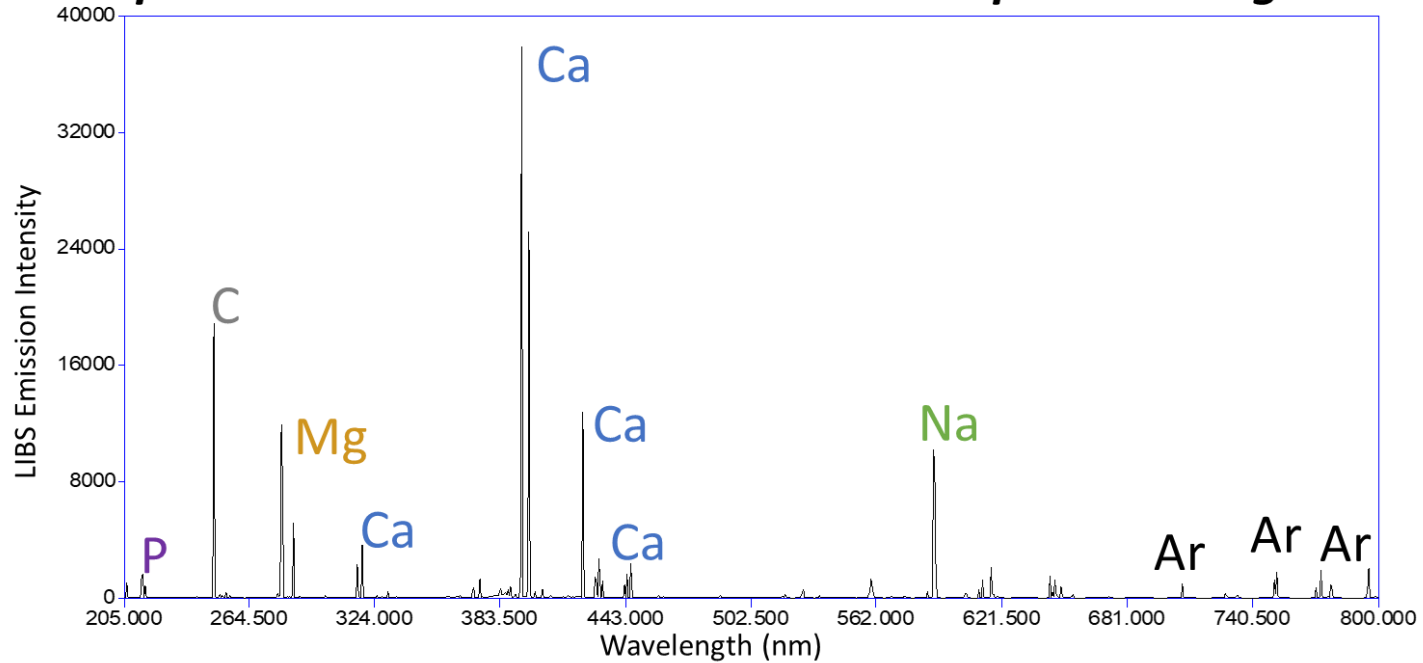


New Apparatus Data

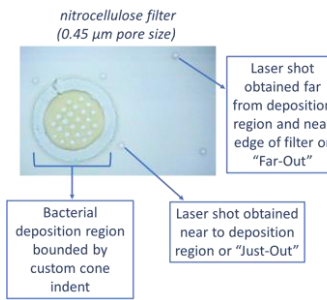

Just-Out and Far-Out LIBS Spectra Overlaid

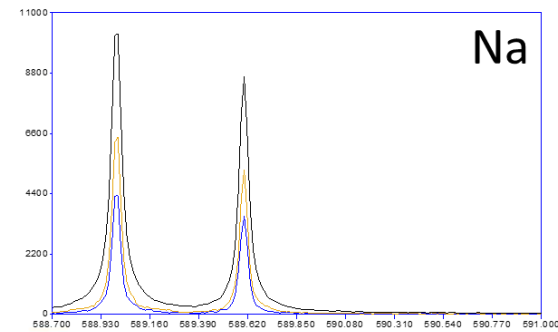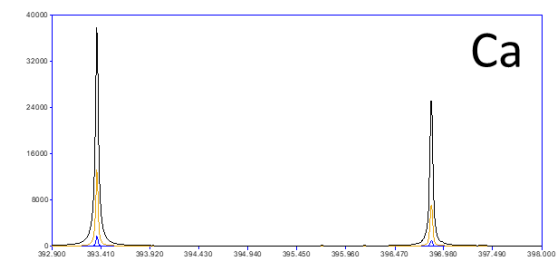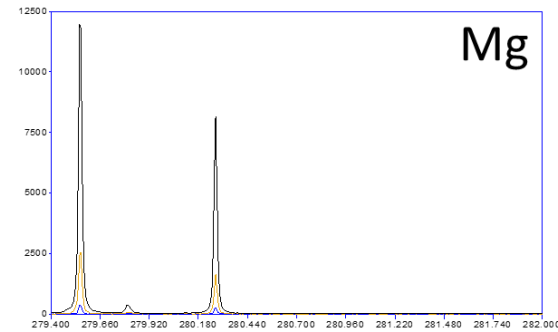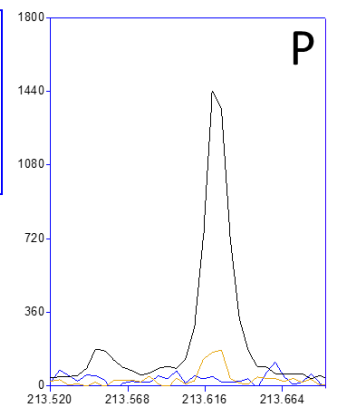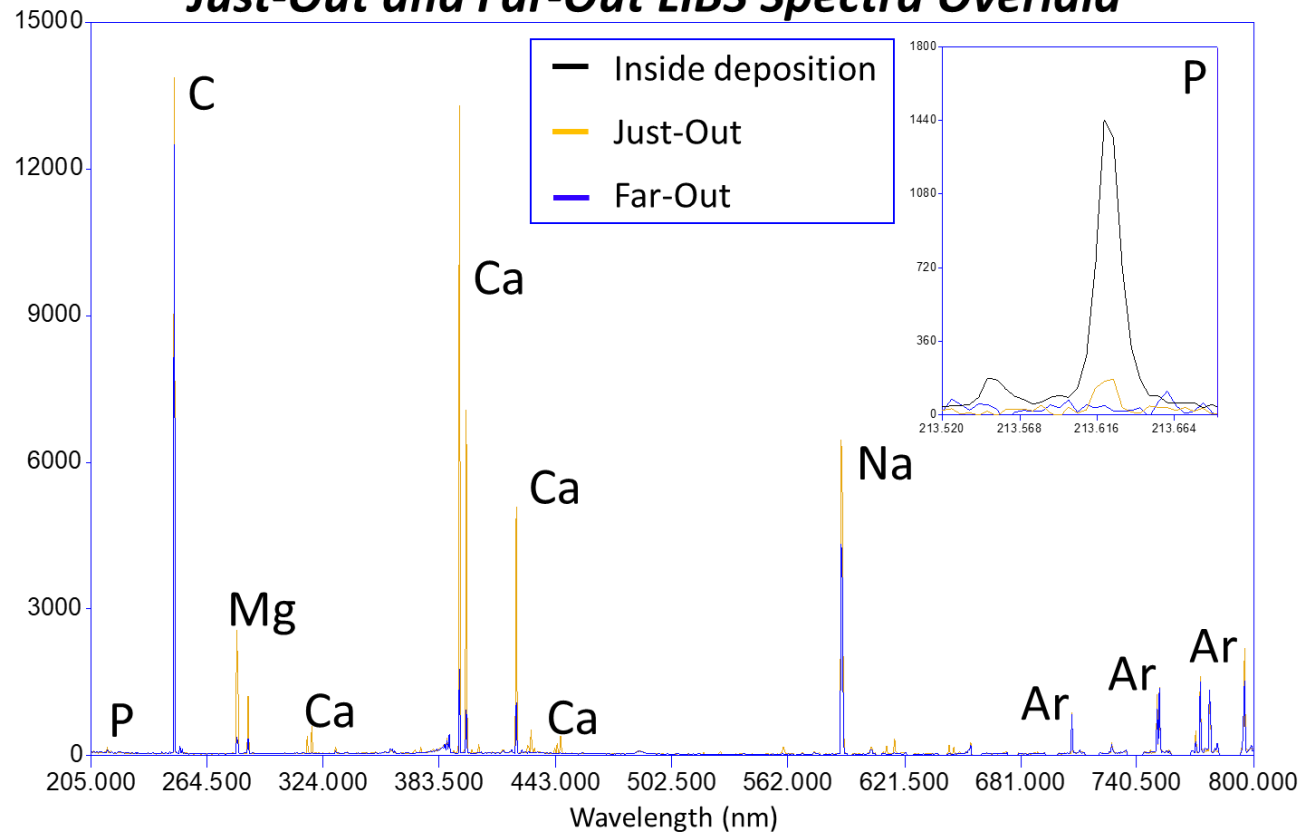

Detecting Infection in aCSF with a PLS-DA

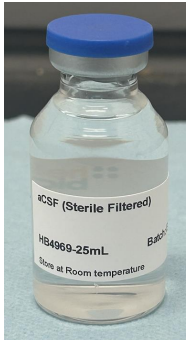

Bacterial-negative aCSF specimens
procured from a biochem supply vendor

Prepared Samples:


- 11 filters each of *S. aureus*,
E. coli, *M. smegmatis*
- 13 filters of sterile aCSF
- 25 LIBS laser shots from each
deposition

DETECTION IN URINE (individual spectra)

99.8% sensitivity – 825 bacterial spectra
100% specificity – 325 urine spectra
11,000 cells per spectrum

DETECTION IN URINE (spectra added)

33 of 33 filters positive, 100% sensitivity
13 of 13 filters negative, 100% specificity

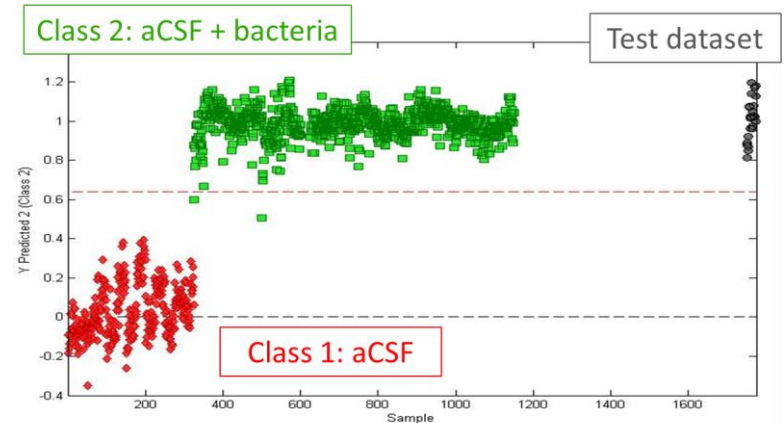

Diagnosing Infection in aCSF with PCA-ANN



Bacterial-negative aCSF specimens procured from a biochem vendor


- Prepared Samples:**
- 11 filters each of *S. aureus*, *E. coli*, *M. smegmatis*
 - 13 filters of sterile aCSF
 - 25 LIBS laser shots from each deposition

Internal cross-validation (80:20 split) randomly selects 20% of the dataset to be tested against remaining 80% of spectra.

80:20 Cross-Validation Test Results				
	<i>E. coli</i>	<i>S. aureus</i>	<i>M. smegmatis</i>	Average
Sensitivity	95.8 %	96.0 %	100.0 %	97.2 %
Specificity	98.2 %	98.0 %	99.7 %	98.6 %
Classification Accuracy	97.4 %	97.3 %	99.8 %	98.2 %

External validation tests individually compare one filter (25 spectra) against left over spectra.

External Validation Test Results				
	<i>E. coli</i>	<i>S. aureus</i>	<i>M. smegmatis</i>	Average
Sensitivity	74.6 %	71.1 %	99.5 %	81.7 %
Specificity	87.6 %	87.9 %	97.0 %	90.9 %
Classification Accuracy	83.3 %	82.3 %	97.9 %	87.8 %

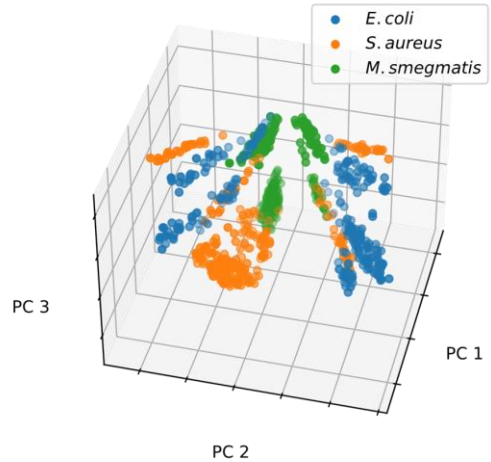

The first three PC scores were plotted to visualize the variance in the data




Conclusions

- Bacteria can be **detected** (PLSDA) and **diagnosed** (PCA-ANN) in blood, urine, and CSF.

PCA-ANN With Full Spectrum Data (Urine)			
	<i>S. aureus</i>	<i>E. coli</i>	<i>E. cloacae</i>
Sensitivity	100 %	100 %	91.67 %
Specificity	100 %	95.83 %	100 %
Classification Error	0.00 %	2.09 %	4.17 %


PCA-ANN With Full Spectrum Data (Blood)				
	<i>S. aureus</i>	<i>E. coli</i>	<i>E. cloacae</i>	<i>P. aeruginosa</i>
Sensitivity	100 %	100 %	100 %	100 %
Specificity	100 %	100 %	100 %	100 %
Classification Error	0.00 %	0.00 %	0.00 %	0.00 %

PCA-ANN With Full Specrum Data (aCSF)				
	<i>E. coli</i>	<i>S. aureus</i>	<i>M. smegmatis</i>	Average
Sensitivity	95.8 %	96.0 %	100.0 %	97.2 %
Specificity	98.2 %	98.0 %	99.7 %	98.6 %
Classification Accuracy	97.4 %	97.3 %	99.8 %	98.2 %

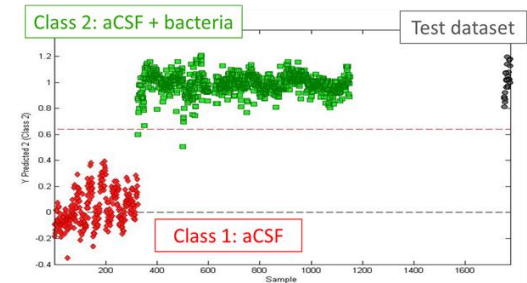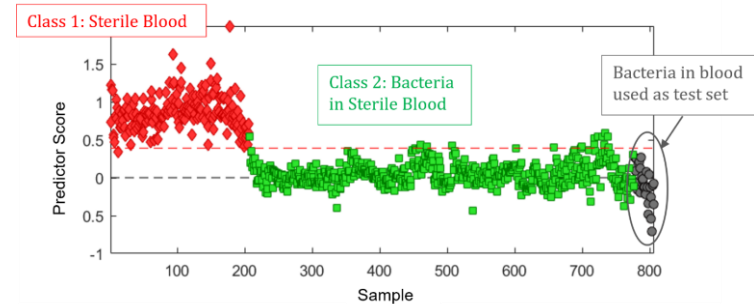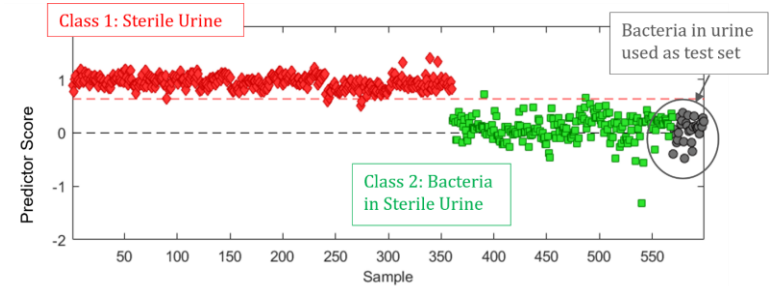


Conclusions

- Adding spectra improves discrimination in PLSDA
 - Detection of bacteria in water, blood, urine, and aCSF improved; **sensitivity = 100%, specificity = 100%**
- Rigorous cleaning of cone & usage of ultrapure water reduces background signal

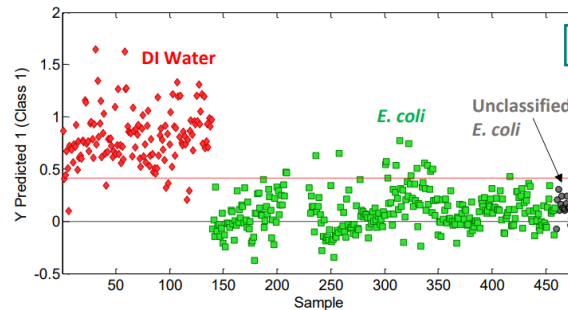

Adding 30 single-shot spectra to create one measurement per filter improves classification in PLSDA

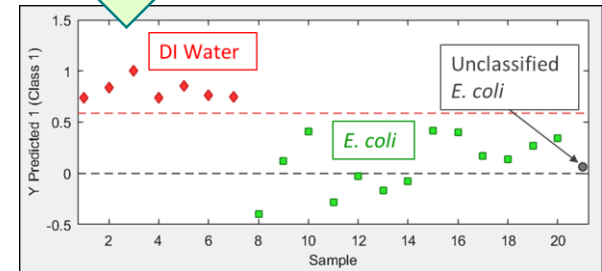

- Using PCA-ANN on full spectrum data provides the best results for discrimination between bacterial species (using 80:20 split)
- Average sensitivity = **94 %**
- Average specificity = **96 %**

	PCA-ANN With Full Spectrum Data		
	E. coli	S. aureus	E. cloacae
Sensitivity	98.04 %	93.27 %	91.23 %
Specificity	97.71 %	97.22 %	96.12 %
Classification Error	2.13 %	4.28 %	6.33 %

Approximate increase from DFA \approx 30 % (sensitivity), 16% (specificity)



Future Work

Improve external validation using PCA-ANN for bacterial species

- New 3D-printed centrifuge insert with integrated concentration cone help to improve signal-to-noise – revisit blood/urine.
- Work to further optimize the PCA-ANN algorithm (adding hidden layers?)
- Determine what is causing the “between-filter” variance in the measurements

Discrimination of lower concentrations of cells to find limit of identification (LOI)

- Engage microbiology students to investigate the behavior of cells in clinical specimens

Extend analysis to other bacterial strains

- *Mycobacterium tuberculosis* and *Streptococcus mitis*



If you want to learn more:

E.J. Blanchette. *Detection and Diagnosis of Bacterial Pathogens in Blood and Urine Using Laser-Induced Breakdown Spectroscopy*. Master's thesis, University of Windsor, 2022.

E.J. Blanchette et al., "Detection and Classification of Bacterial Cells After Centrifugation and Filtration of Liquid Specimens Using Laser-Induced Breakdown Spectroscopy," *Applied Spectroscopy* **76**, 2022, pp. 894-904. <https://doi.org/10.1177/000370282210927>

E.J. Blanchette, E.A. Tracey, A. Baughan, G.E. Johnson, H. Malik, C.N. Alionte et al. "Detection and diagnosis of bacterial pathogens in blood using laser-induced breakdown spectroscopy," *Spectrochim. Acta B* **215** (2024) 106911. <https://doi.org/10.1016/j.sab.2024.106911>

E.J. Blanchette, E.A. Tracey, A. Baughan, G.E. Johnson, H. Malik, C.N. Alionte, et al. "Detection and diagnosis of bacterial pathogens in urine using laser-induced breakdown spectroscopy," *Spectrochim. Acta B* **216** (2024) 106944. <https://doi.org/10.1016/j.sab.2024.106944>





Dr. Steven Rehse

- Home
- Latest Developments
- Students
- Research Activities
- Publications**
- Conferences & Presentations
- DoD STEM Virtual Great Lakes Region and Tribal Nations Sessions
- Archive
- About Me / Bio

UWindsor Home > Dr. Steven Rehse > Publications


Publications

- > Book Chapters
- > **Journal Articles**
- > Dissertations




Acknowledgements for the people who did the work...

Emma

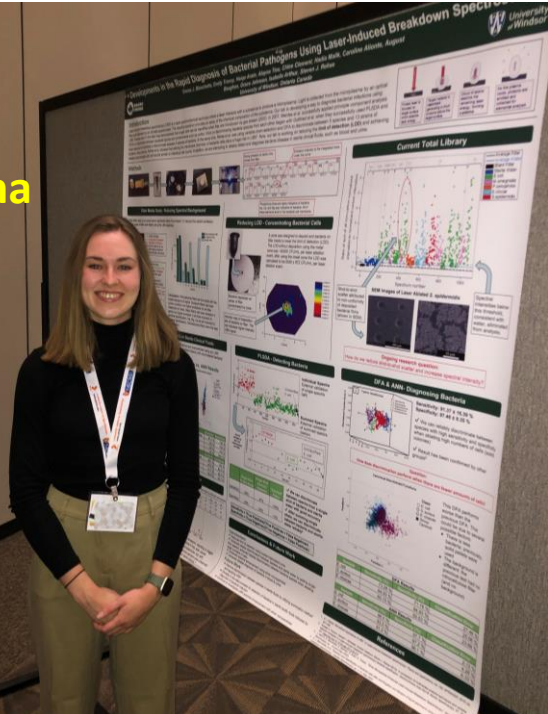


Grace

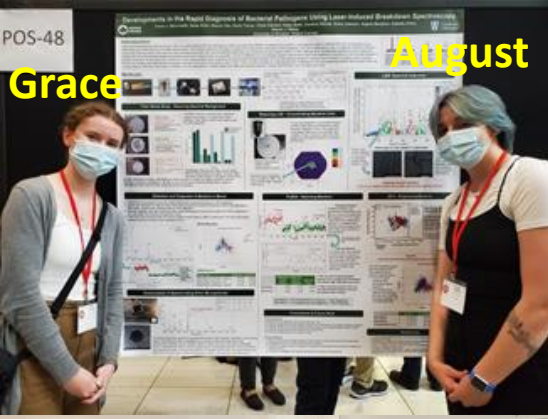

August


Emily

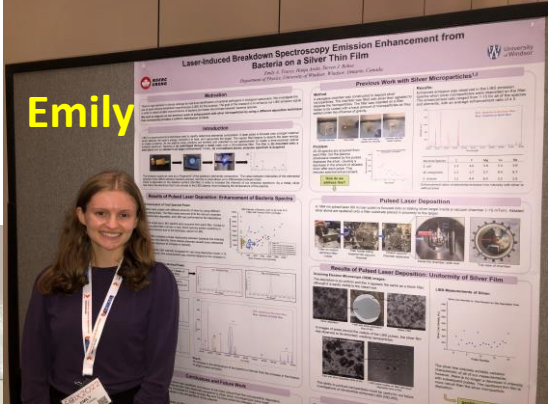


Abdullah

Caroline


Nick

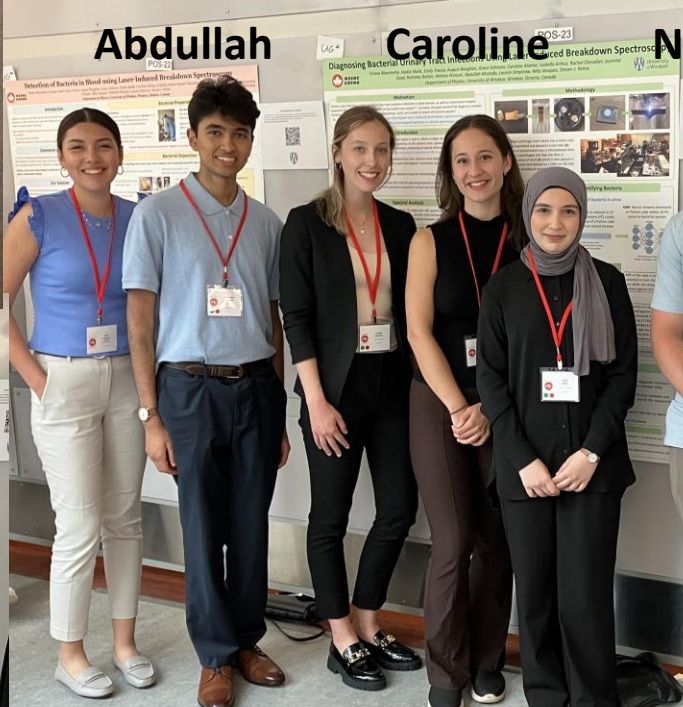

Izzy

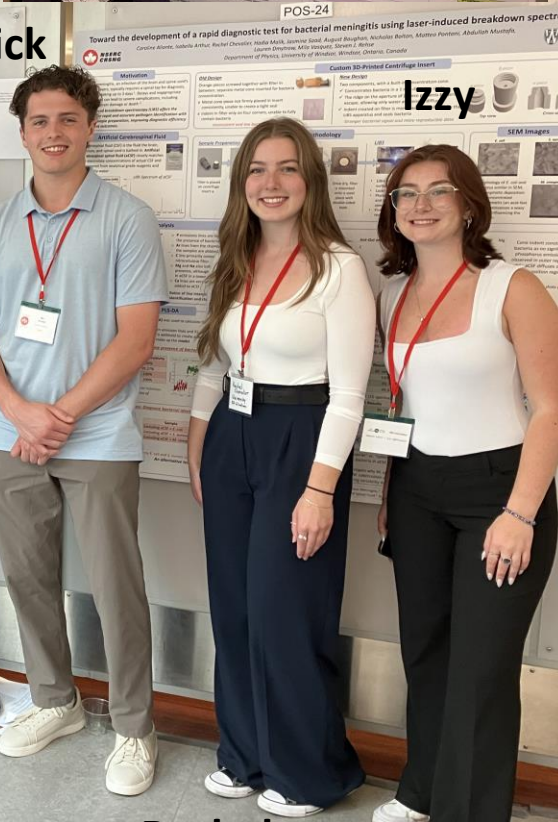

Jasmine

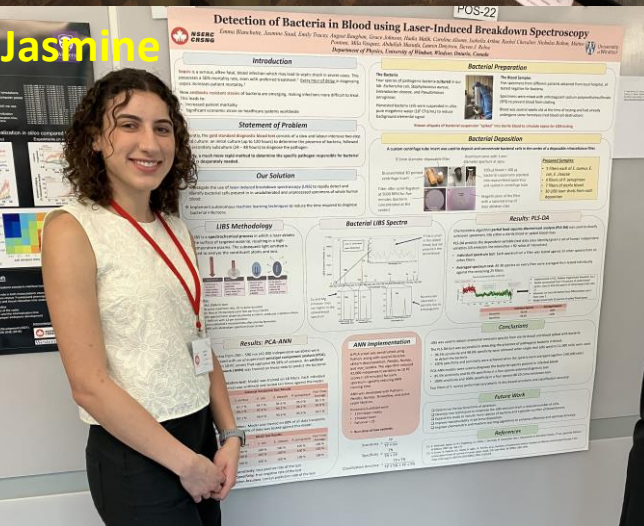

Mila

Lauren

Hadia

Rachel

Funding and Acknowledgements

We gratefully acknowledge funding for this project provided by:

- A **Natural Sciences and Engineering Research Council of Canada** Discovery grant and a Research Tools and Instruments grant
- A **Canada Foundation for Innovation** Leaders Opportunity Fund grant
- An **Ontario Research Fund** Small Infrastructure Funds grant
- **University of Windsor** Outstanding Scholars program
- **University of Windsor** Faculty of Science



Thank you!



YouTube




Rehse UWindsor Class Channel

@rehseuwindsorclasschannel6484 222 subscribers 52 videos




More about this channel >

Webpage: www.uwindsor.ca/people/rehse/

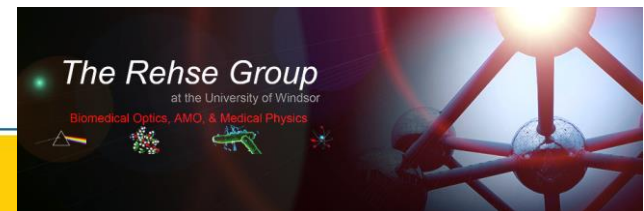





2) removal of samples mass (ablation)

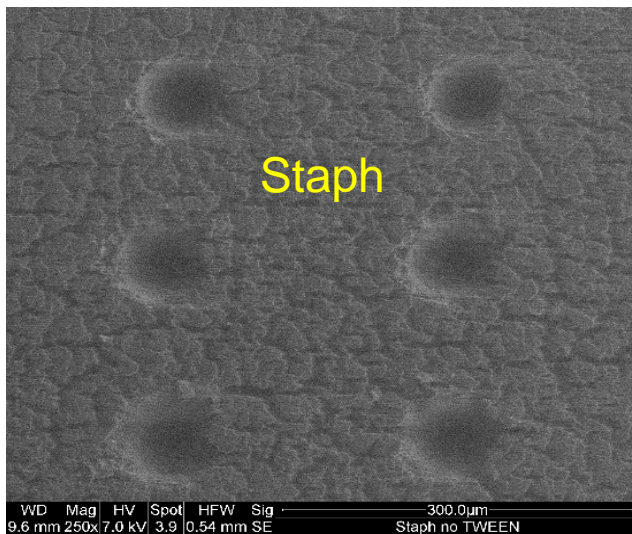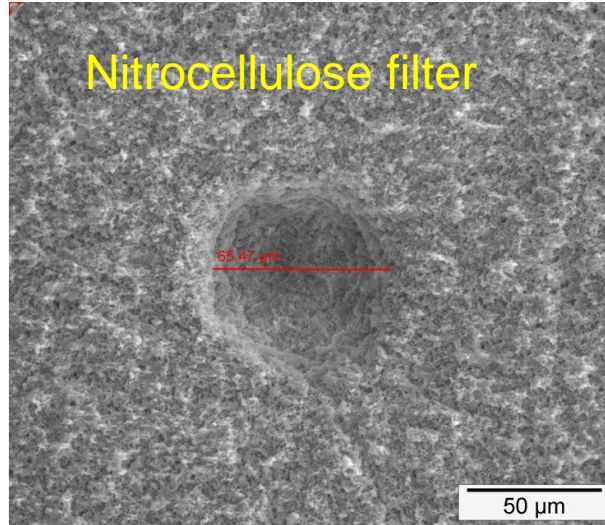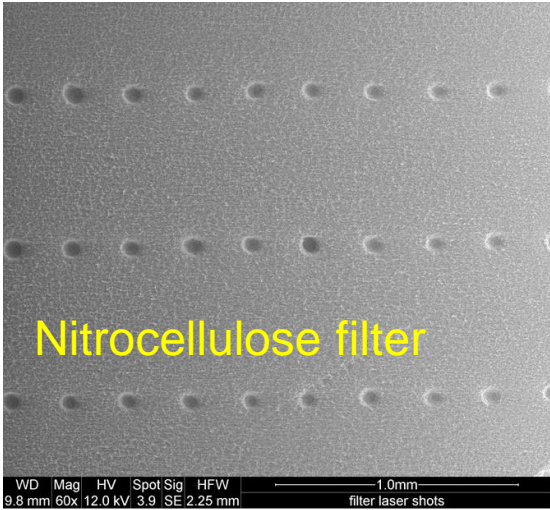



M. smegmatis

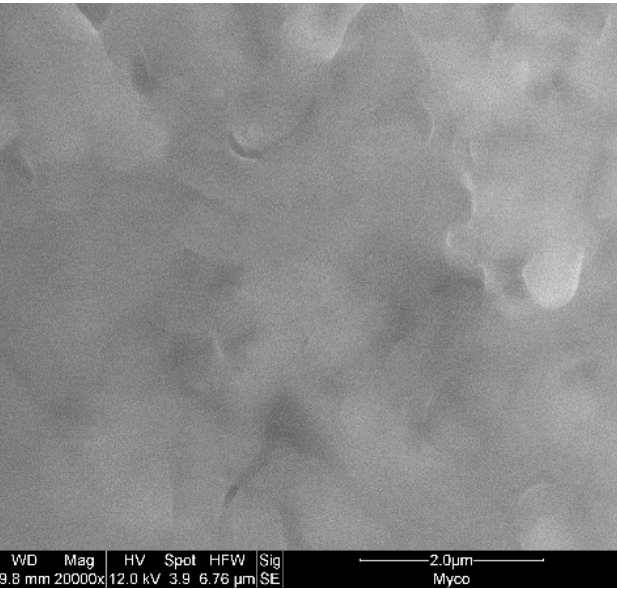

E. coli

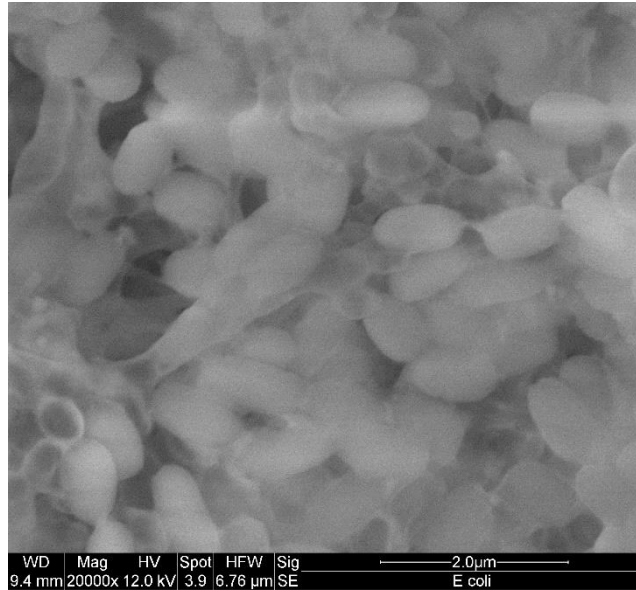


S. aureus

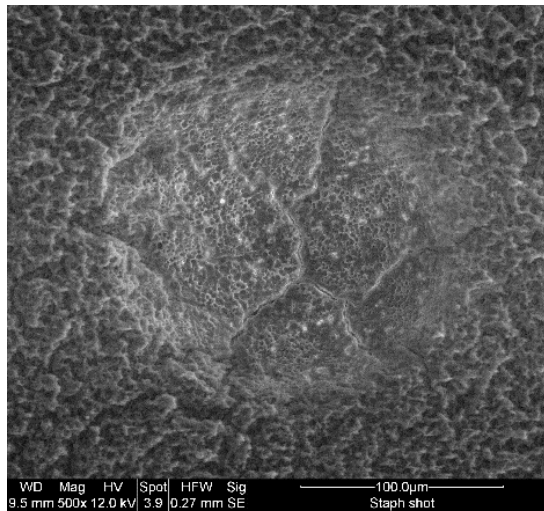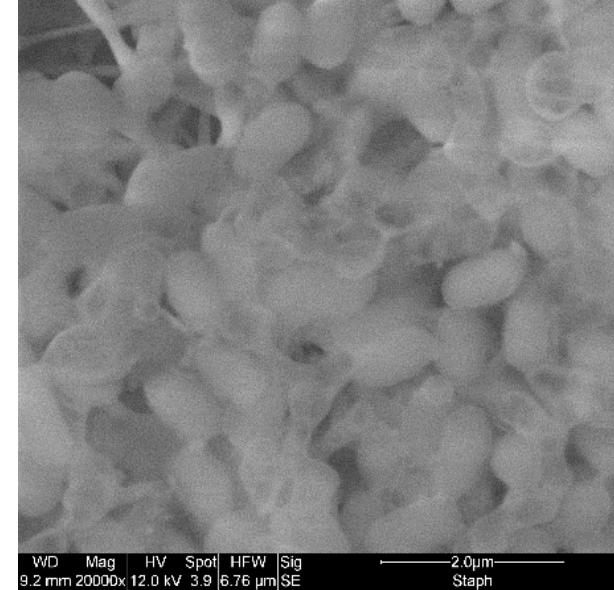

PCA-ANN Results – Identifying Bacteria


GOAL: Classify the species of bacteria in urine

Input: Entire LIBS spectrum from 200 nm to 590 nm (42,000 variables)

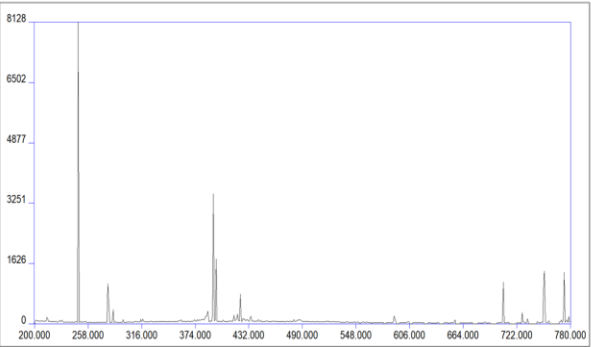


12 spiked urine filters x 30 spectra/filter

PCA: Spectrum is reduced to 10 principal components (PC) scores, through the use of a Python code developed from several libraries

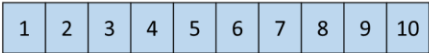

Python libraries used: *sklearn.decomposition, Pandas, Numpy, and mpl_toolkits*

ANN: Neural network developed on Python code relates 10 PC scores to bacterial species

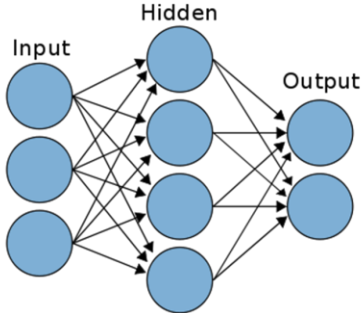


Python libraries used: *Pandas, Numpy, Tensorflow, and Scikit-Learn*





**4) expansion and
element specific
emission (atomic or
ionic)**

**spontaneous emission
as atoms/ions decay to
ground state**

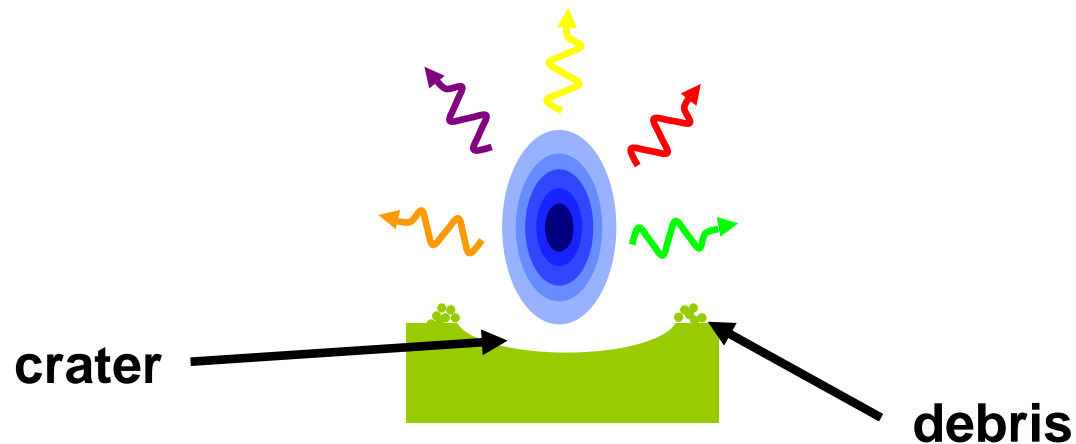


- The dynamic evolution of the plasma plume is then characterized by a fast expansion and subsequent cooling.
- Approximately 1 microsecond after the ablation pulse, spectroscopically narrow atomic/ionic emissions may be identified in the spectrum.

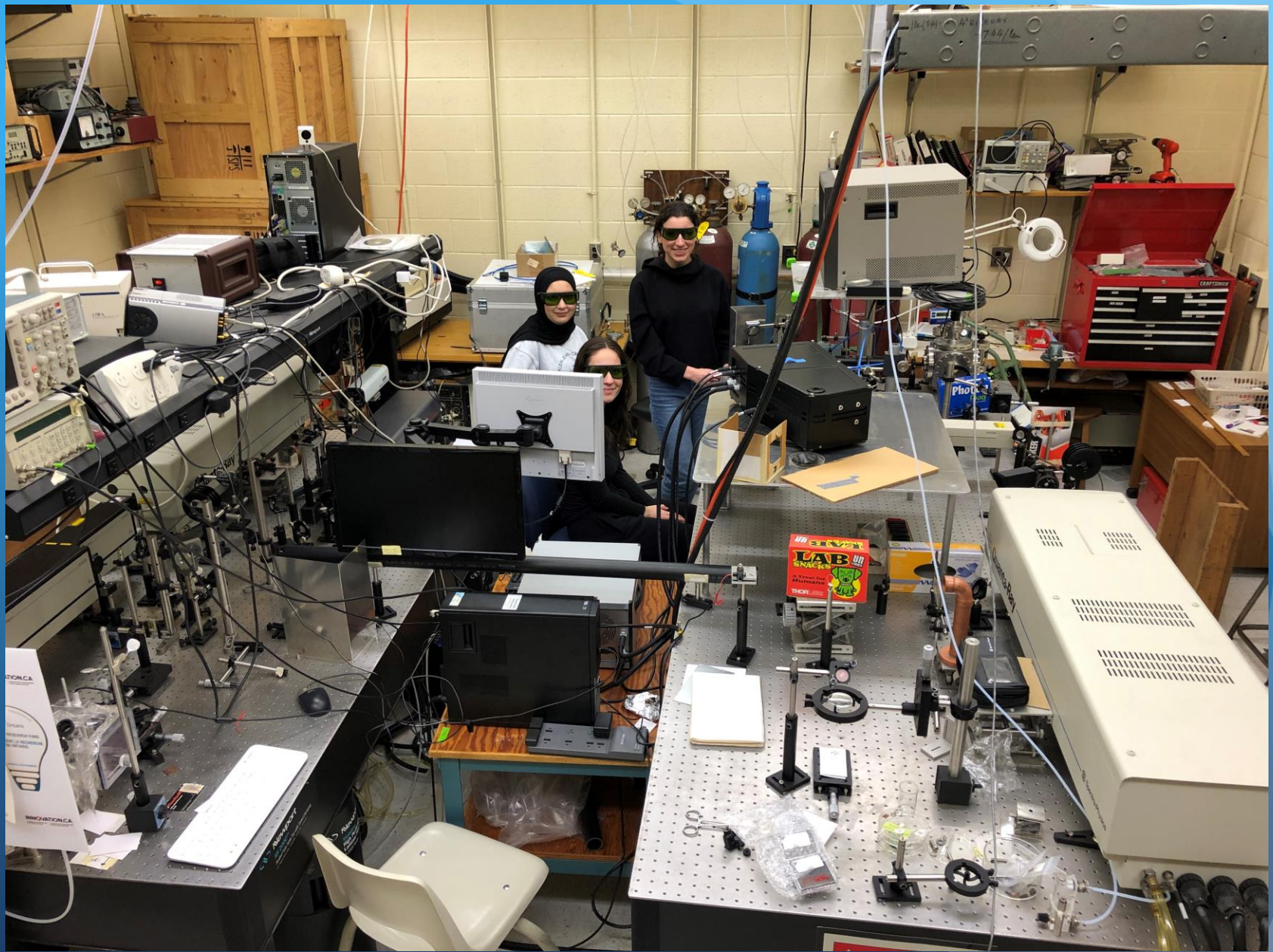

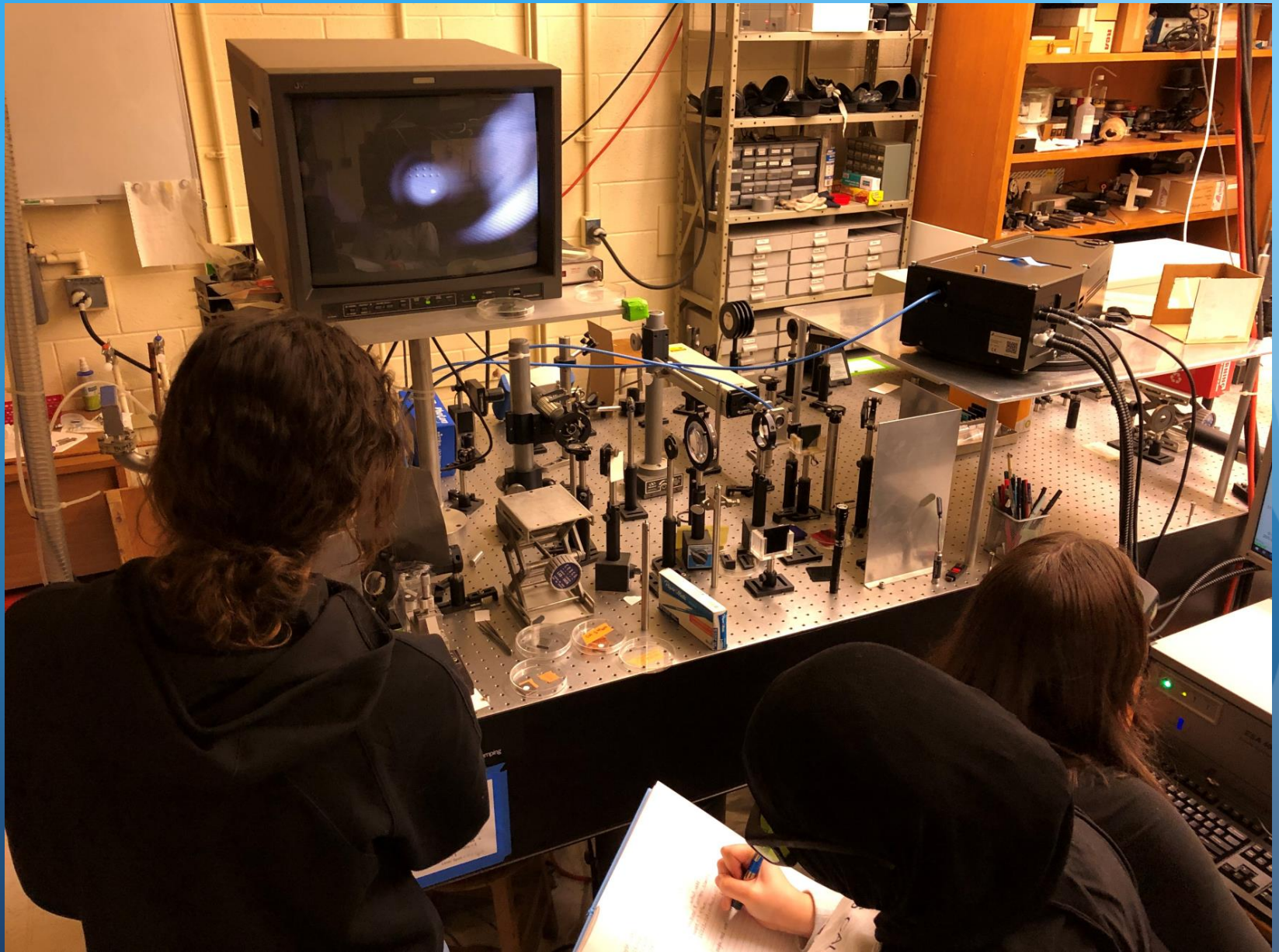


The Centrifuge Insert

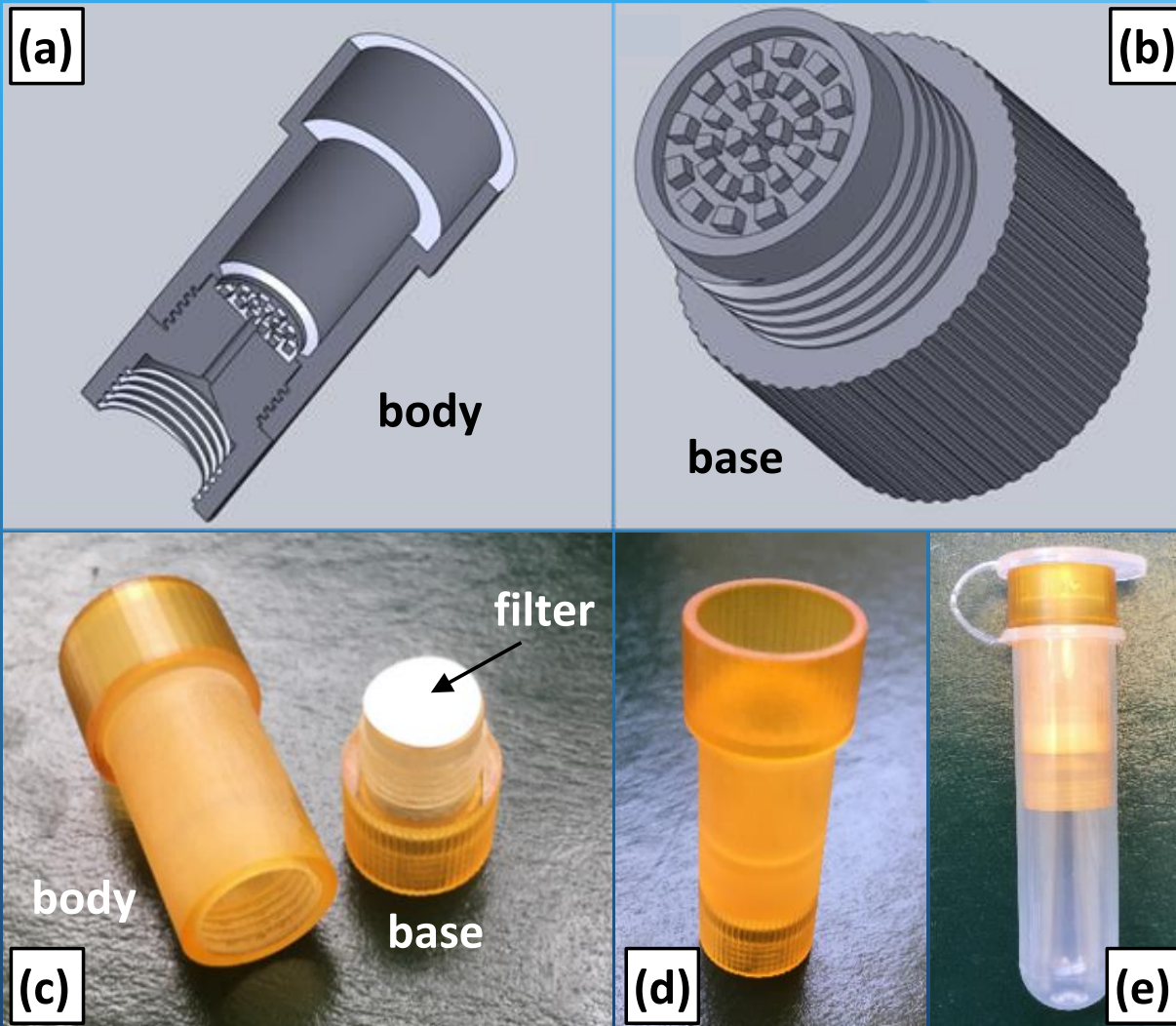

Concentrating Bacteria With a Cone

19 mm long Al cone

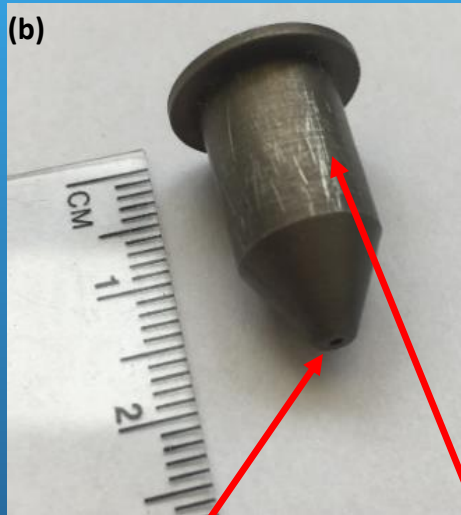

Centrifuge tube cap presses cone into filter



1 mm hole at apex




Holds 1 mL of fluid




Cone vertex press fit into filter

Shooting Bacteria Concentrated With Cone

Total LIBS intensity (A.U.)

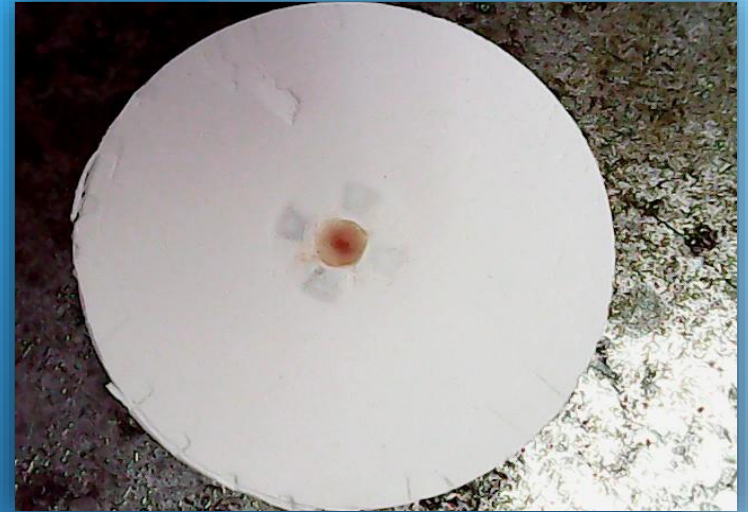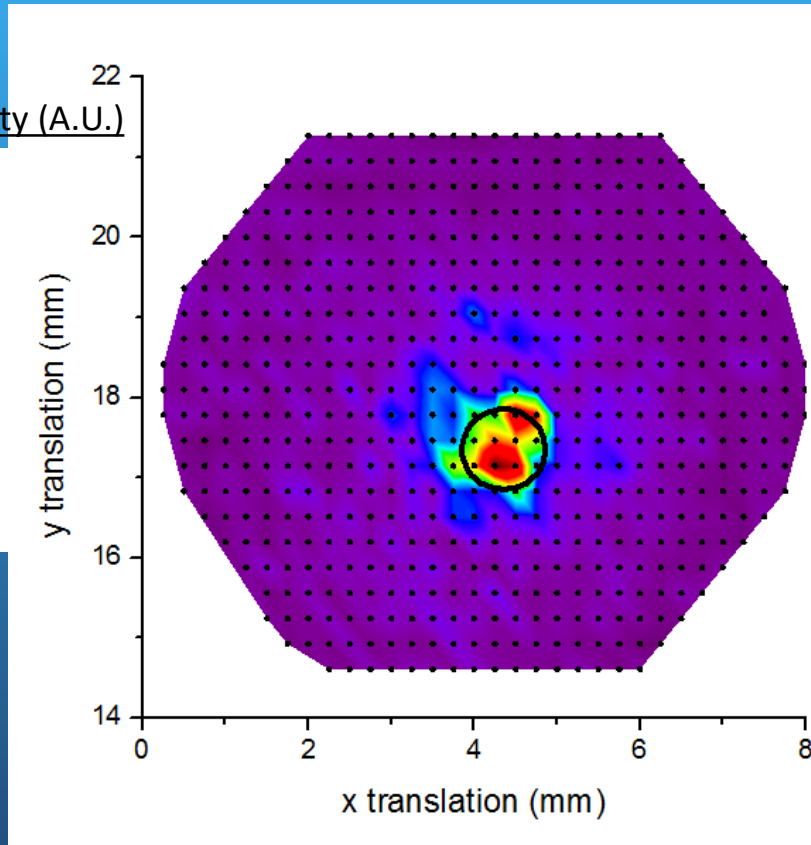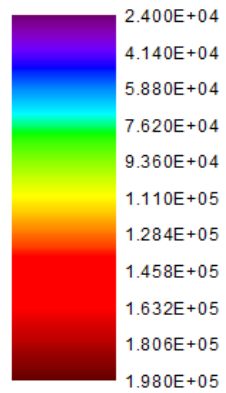

Limits of Detection

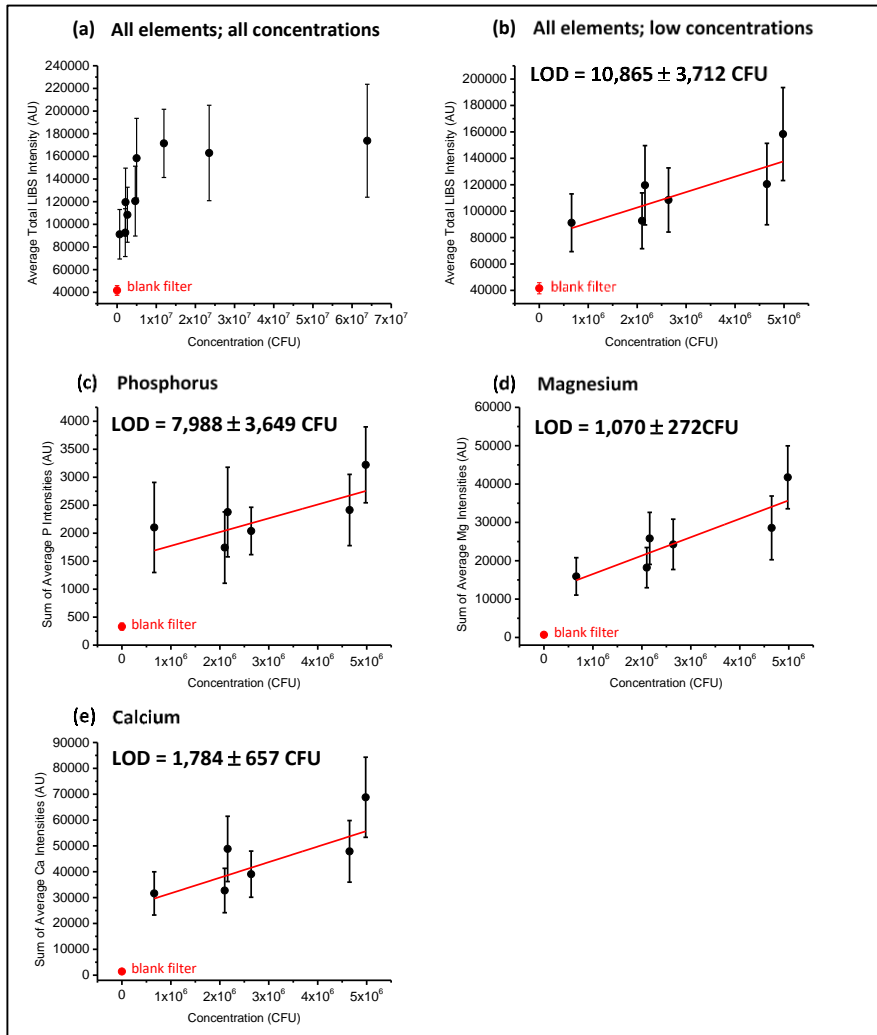

A calibration curves constructed from forty spectra obtained from each of nine different concentrations.

LIBS bacterial limit of detection of $10,865 \pm 3,712$ CFU per laser ablation event for bacteria deposited on filters using the metal cone.



LOD's calculated using only certain elements observed in the LIBS spectra and present in very low concentrations in the filter were even lower:

$1,070 \pm 272$ CFU for magnesium
 $1,784 \pm 657$ CFU for calcium.

LOD on filter better, but number of cells required in fluid specimen is **WAY** lower!

When performed with no background

- ✓ We can identify a bacterial species, certainly its genus, with high sensitivity and specificity (confirmed by others).
- ✓ We can differentiate strains of *E. coli* (demonstrated by others in MRSA).
- ✓ Multiple multivariate techniques effective at discriminating spectra.

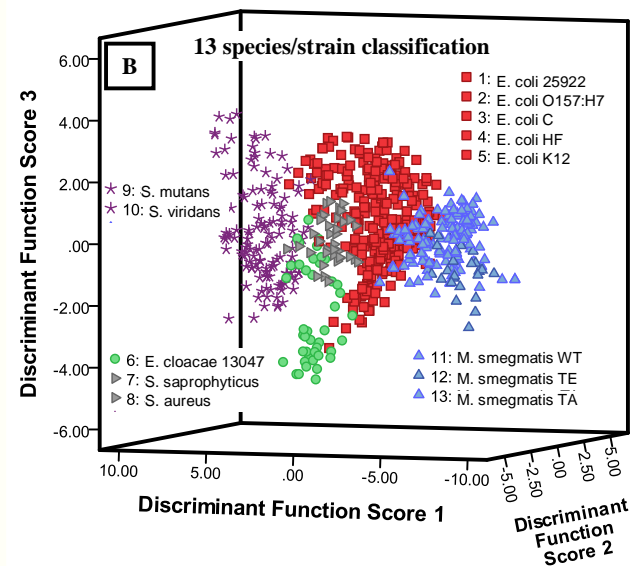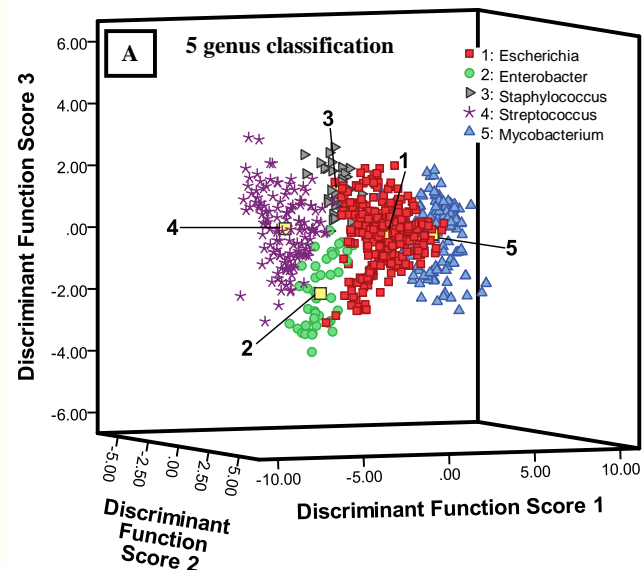

PLSDA

<i>E. COLI</i>	True	False
Positive	95.65%	9.17%
Negative	90.83%	4.35%
<i>STAPHYLOCOCCUS</i>	True	False
Positive	54.05%	0.51%
Negative	99.49%	45.95%
<i>STREPTOCOCCUS</i>	True	False
Positive	95.59%	1.02%
Negative	98.98%	4.41%
<i>MYCOBACTERIUM</i>	True	False
Positive	88.31%	1.06%
Negative	98.94%	11.69%

DFA

<i>E. COLI</i>	True	False
Positive	89.63%	15.95%
Negative	84.05%	10.37%
<i>STAPHYLOCOCCUS</i>	True	False
Positive	86.49%	5.85%
Negative	94.15%	13.51%
<i>STREPTOCOCCUS</i>	True	False
Positive	99.26%	13.32%
Negative	88.68%	0.74%
<i>MYCOBACTERIUM</i>	True	False
Positive	96.10%	4.08%
Negative	95.92%	3.90%

DFA: Sensitivity: $91.37 \pm 16.39 \%$ Specificity: $97.46 \pm 9.35 \%$
PLSDA: Sensitivity: $93.13 \pm 10.25 \%$ Specificity: $90.60 \pm 21.33 \%$

Results: We have already demonstrated...

- LIBS spectral fingerprint is a *sensitive* and *specific* (high rates of true positives, low rates of false positives) test to identify an unknown bacterial specimen or to differentiate between possible identifications
- This spectral fingerprint is *robust* and *reliable*, and exists through time (multiple tests spanning years on same strains of bacteria)
- In addition...



Results: We have already demonstrated...

LIBS spectral fingerprint is:

- growth-medium independent
- independent of state of growth (how “old” the bacteria are)
- independent of whether the bacteria are live or dead (or inactivated by UV light)
- obtainable even when other types of bacteria or contaminants are present (mixed samples)
- obtainable from urine specimens
- capable of strain discrimination
- obtainable from about 500 bacteria

9 publications in Applied Physics Letters, Journal of Applied Physics, Applied Optics, Applied Spectroscopy, Spectrochimica Acta B, and others – confirmed by multiple other groups

(Poster 344) **Bacterial Mounting and Concentration Techniques to Translate Laser-Induced Breakdown Spectroscopy into a Clinical Setting;** Alexandra Paulick

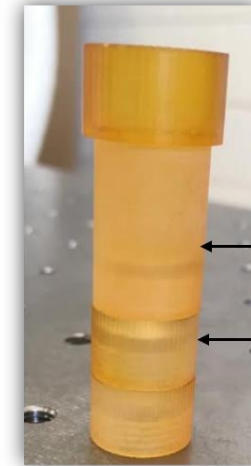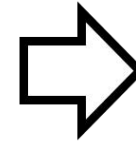

5 μm filter



0.45 μm filter